



NEW!

Stretcher Series

7500

transport / trauma stretcher



PRODUCT CATEGORIES

STRETCHERS

- 7500 Stretcher
- 5400 Universal Stretchers
- 5100 General Stretchers
- Pediatric Stretcher Cribs
- 170 General Stretchers

CARTS

- Surgical Case Carts
- Utility Carts
- Supply Carts
- Distribution Carts
- Multi-Purpose Carts
- Accessories

IV STANDS

- Stainless Steel and Chrome
- Infusion Pump Stand

NURSERY

- Labor & Delivery
- Bassinets

STANDS/KICKBUCKETS

- Mayo Stands
- Basin Stands & Kick Buckets
- Linen Hampers

STOOLS/FOOTSTOOLS

- Seating
- Footstools

TABLES

- Instrument Assembly Station
- The Space Station
- Instrument Tables
- Utility Tables
- Central Supply Work Tables
- Specialty Operating Room Tables

WARMING CABINETS

- Blanket & Fluid Warmers
- Combination Units

MISCELLANEOUS

- Operating Room Cabinets
- Anesthesia and Treatment Cabinets
- Physician Office Equipment
- Sterile Processing Equipment

Elite Series

Blanket Warming Cabinets



P-2220



P-2205



800.822.3501 • www.pedigo-usa.com

4000 SE Columbia Way • Vancouver, WA 98661

09/12 MM-02



operating room cabinets



P-8230 Free-standing cabinet

FEATURES

- Stainless steel cabinetry to maximize storage and space utilization
- Available as free-standing, built-in, or caster base models
- Double panel doors with FDA approved polyethylene foam for glass shock resistance
- $\frac{3}{16}$ " tempered glass doors with stainless steel lockable paddle latch
- Radius glass panels for ease of cleaning
- Adjustable stainless steel door hinges
- Adjustable nylon floor glides for leveling
- Stainless steel shelves adjustable in $\frac{1}{2}$ " increments
- Adaptable to custom installations
- Keyed door latch



Built-in cabinets with trim bezel

OPTIONS

- Stainless steel or Polyvinyl pegboard
- Built-in designs
- Trim bezel for built-in applications
- Caster base with locking casters
- Solid doors
- Additional shelves



Sloped top with caster base

MANY SIZES AND STYLES TO CHOOSE FROM:

Part No.	Dimensions	Shelves	Top
P-8010	25 $\frac{3}{4}$ "W x 18"D x 60"H	3	Flat
P-8015	25 $\frac{3}{4}$ "W x 18"D x 60"H	3	Sloped
P-8020	25 $\frac{3}{4}$ "W x 24"D x 60"H	3	Flat
P-8025	25 $\frac{3}{4}$ "W x 24"D x 60"H	3	Sloped
P-8030	47"W x 18"D x 60"H	3	Flat
P-8035	47"W x 18"D x 60"H	3	Sloped
P-8040	47"W x 24"D x 60"H	3	Flat
P-8045	47"W x 24"D x 60"H	3	Sloped
P-8110	25 $\frac{3}{4}$ "W x 18"D x 72"H	4	Flat
P-8115	25 $\frac{3}{4}$ "W x 18"D x 72"H	4	Sloped
P-8120	25 $\frac{3}{4}$ "W x 24"D x 72"H	4	Flat
P-8125	25 $\frac{3}{4}$ "W x 24"D x 72"H	4	Sloped

Part No.	Dimensions	Shelves	Top
P-8130	47"W x 18"D x 72"H	4	Flat
P-8135	47"W x 18"D x 72"H	4	Sloped
P-8140	47"W x 24"D x 72"H	4	Flat
P-8145	47"W x 24"D x 72"H	4	Sloped
P-8210	25 $\frac{3}{4}$ "W x 18"D x 84"H	5	Flat
P-8215	25 $\frac{3}{4}$ "W x 18"D x 84"H	5	Sloped
P-8220	25 $\frac{3}{4}$ "W x 24"D x 84"H	5	Flat
P-8225	25 $\frac{3}{4}$ "W x 24"D x 84"H	5	Sloped
P-8230	47"W x 18"D x 84"H	5	Flat
P-8235	47"W x 18"D x 84"H	5	Sloped
P-8240	47"W x 24"D x 84"H	5	Flat
P-8245	47"W x 24"D x 84"H	5	Sloped





ANESTHESIA AND
TREATMENT CABINETS



P-1100-SS

ANESTHETIST TABLE WITH DRAWER

OVERALL SIZE: 21³/₈" W x 14" D x 34¹/₄" H

TRAYS: 12⁵/₈" W x 5⁵/₈" D x 19⁵/₈" H

DRAWER: 15³/₈" W x 12¹/₂" D x 3¹/₄" H

P-1101-SS

Same as P-1100-SS without drawer

FEATURES

- Ideal cart for storage of instruments, gloves, sutures, etc.
- Two removable stainless steel trays
- Towel bar and two mask hooks
- All welded and polished stainless steel construction
- 2" double ball bearing swivel casters



P-1105-SS

ANESTHETIST CABINET WITH DRAWERS

OVERALL SIZE: 19¹/₈" W x 15⁵/₈" D x 33³/₈" H

UPPER 3 DRAWERS: 13¹/₂" W x 12" D x 3¹/₂" H

BOTTOM DRAWER: 13¹/₂" W x 12" D x 7" H

FEATURES

- Versatile anesthetist cabinet provides extra storage space
- Full extending drawer hangers
- Removable 12⁵/₈" x 19¹/₈" stainless steel tray
- All welded & polished stainless steel construction
- 3" double ball bearing swivel casters



P-5090

SMALL TREATMENT CABINET

OVERALL SIZE: 16" W x 20" D x 33" H

DRAWER: 13¹/₂" W x 12" D x 3¹/₂" H

FEATURES

- Solid stainless steel retaining rail
- Internal shelf for optimal storage
- Optional locking rod
- All welded and polished stainless steel construction
- 2" double ball bearing swivel casters



P-5091

LARGE TREATMENT CABINET

OVERALL SIZE: 28" W x 18" D x 33" H

TOP DRAWER: 21¹/₄" W x 16" D x 3¹/₄" H

BOTTOM DRAWER: 21¹/₄" W x 16" D x 5³/₈" H

FEATURES

- Solid stainless steel retaining rail
- Internal shelf for optimal storage
- Optional locking rod
- All welded and polished stainless steel construction
- 2" double ball bearing swivel casters



P-7202

CARDIAC AND ANESTHETIST CABINET

OVERALL SIZE: 28" W x 18" D x 33" H

TOP DRAWER: 21¹/₄" W x 16" D x 3¹/₄" H

BOTTOM DRAWER: 21¹/₄" W x 16" D x 5³/₈" H

P-7202-C

Same as P-7202 with complete accessory package


OPTIONAL ACCESSORIES FOR P-7202

- O2 holder (holds 2 E-size tanks)
- I.V. rod
- Electrical outlet with 25' cord
- Push handle
- Cabinet locking rod

FEATURES

- Designed for cardiac and other emergency uses or as a large anesthetist cabinet
- Provides space for instruments, vials, and other necessary equipment
- All welded and polished stainless steel construction
- 4" double ball bearing swivel casters

shown with optional accessories

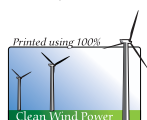



SURGICAL CASE CARTS



CDS-233 shown with optional shelves

C D S - 2 3 3

CABINET SIZE

23"W x 24¼"D x 47"H

OVERALL SIZE

28⁵/₈"W x 28⁵/₈"D x 56"H

C D S - 2 3 5

CABINET SIZE

28½"W x 28½"D x 47"H

OVERALL SIZE

34¹/₈"W x 32³/₈"D x 56"H



CDS-242 shown with optional shelf

C D S - 2 4 2

CABINET SIZE

40¼"W x 24¼"D x 30"H

OVERALL SIZE

45⁷/₈"W x 28⁵/₈"D x 39"H



CDS-245 shown with optional shelves

C D S - 2 4 5

CABINET SIZE

40¼"W x 24¼"D x 47"H

OVERALL SIZE

45⁷/₈"W x 28⁵/₈"D x 56"H



CDS-256-MS with optional leaf shelf

C D S - 2 5 6 - M S

CABINET SIZE

23"W x 24¼"D x 30"H

OVERALL SIZE

28⁵/₈"W x 28⁵/₈"D x 39"H

C D S - 2 5 6 - L D

CABINET SIZE

23"W x 24¼"D x 21"H

OVERALL SIZE

28⁵/₈"W x 28⁵/₈"D x 30"H

FEATURES

- All carts built to Steril-Gard specifications
- Fully-welded all stainless steel construction
- Sound-proofed tops and doors for quiet operation
- Exterior push handles for better control and more interior space
- One-piece non-marring continuous wrap-around bumper
- Positive lock door latches
- Sealed doors for sterility
- All carts can be manufactured to individual specifications
- Quick delivery

OPTIONS

STAINLESS STEEL CASTERS

- Ideal for carts that are cart washed daily
- Large 6" easy rolling casters
- All stainless steel material and bearings for a rust-free caster
- Two lubrication points for ease of maintenance

WIRE SHELVES

- Electro-polished stainless steel
- Easy to adjust or remove
- ³/₁₆" diameter wire welded to ¼" diameter wire support rods
- Maximum opening between rods: 1" x 6½"
- Option of roll-out or stationary shelves

SOLID SHELVES

- Steril-Gard construction
- Easy to adjust or remove
- 18-gauge stainless steel
- Option of roll-out or stationary shelves

For a complete list of optional accessories, please view our "Closed Cart Accessories" brochure.





Pedigo is proud to introduce our sealed Surgical Case Carts. These carts provide an added level of protection against dust and debris. Features a full seal around the door perimeter and center stop. Withstands cart wash/dry temperatures and sanitizing detergents.

OPTIONAL FEATURES

- Security Latch
- Heavy Duty 6" x 2" Casters
- Color Identification Wheel
- Sloped Bottom
- Barcode Plate Holder
- 8½" x 11" Card Holder

Color Identification Wheel

As an added option, a color ID wheel allows the hospital to identify various case carts for specific applications. This allows for easy identification of completed carts without breaking the tamper proof seal and compromising the contents.



Sloped Bottom

The 2° slope forces fluids to the back of the cart for fluid containment. The contained fluids will drain out the front of the cart during cart washing.

Heavy Duty Casters

Optional Stainless Steel 6" diameter casters with a 2" wide wheel. Higher load rating for improved durability. Two Swivel Casters on the pushing end of the cart and two rigid casters on the opposite end.



Standard Casters

Carts come with 6" diameter swivel chrome casters and a 1¼" wide wheel.



Sealed Doors

Integrated silicon door seal installed on the interior of the doors, including the seam between the doors. The seal protects the sterile contents inside of the cart from the outside environment.



Locking Shelves

Optional lock on roll out shelf holds the shelf in either the open or closed position.



Total Lock Brake

Optional feature recommended for offsite sterile processing. Casters can be locked after carts are loaded onto trucks.





STERILE PROCESSING
EQUIPMENT



shown with hospital accessories

C D S - 1 6 0

INSTRUMENT CONTAINER WASH CART

OVERALL SIZE

49"W x 26"D x 67"H

FEATURES

- Stainless steel construction
- Fully cart-washable
- Built to Steril-Gard specifications
- Fully welded frame
- Donut bumpers on legs
- 5" diameter stainless steel swivel casters
- Five stainless steel wire shelves
- Six adjustable container lid dividers



shown with hospital accessories

C D S - 1 6 2

STERILE PROCESSING WASH CART

OVERALL SIZE

49"W x 26"D x 66½"H

FEATURES

- Stainless steel construction
- Fully cart-washable
- Built to Steril-Gard specifications
- Fully welded frame
- Donut bumpers on legs
- 5" diameter stainless steel swivel casters
- Three stainless steel wire shelves with retaining rods
- Holds instrument containers, tote boxes, and basins



shown with optional accessories

C D S - 3 2 6 0

CENTRAL SUPPLY WORK TABLE

OVERALL SIZE

60"W x 32"D x 34"H

FEATURES

- Stainless steel construction
- Height adjustable floor glides

OPTIONAL ACCESSORIES

- 4" diameter swivel casters with brakes
- Sterile Wrap hangers measure 56"W x 30½"H
- Full-length lower shelf
- Tote box drawers, 5" or 10" deep; or stainless steel drawers



Caster with brake





SUPPLY CARTS



CDS-152-12

TOTE BOX CART

OVERALL SIZE

36¼"W x 24⅛"D x 45"H

FEATURES

- Ideal cart for storage of instruments, gloves, sutures, etc.
- All welded stainless steel construction, single piece frame
- Outboard stem-mounted 5" casters for stability
- Wrap around bumper mounted away from the frame allows easy cleaning
- Includes 12 standard 5" deep tote boxes or 6 large 10" tote boxes or any combination of the two sizes



CDS-152-24

TOTE BOX CART

OVERALL SIZE

51½"W x 24⅛"D x 57"H

FEATURES

- All welded stainless steel construction, single piece frame
- Outboard stem-mounted 5" casters with diagonally mounted swivel locks provide stability
- Wrap around bumper mounted away from frame for easy cleaning
- Includes 24 standard 5" deep tote boxes or 12 large 10" deep tote boxes or any combination of the two sizes



shown with optional accessories

CDS-177

CENTRAL SUPPLY CART

OVERALL SIZE

55⅜"W x 26⅛"D x 58"H

FEATURES

- All welded stainless steel construction for one single piece frame
- Outboard plate-mounted 8" casters with diagonally mounted swivel locks provide stability
- Wrap around bumper mounted away from frame for easy cleaning
- Includes seven standard 5" deep tote boxes or up to three large 10" deep tote boxes
- Optional shelves and accessories available to configure the cart to suit individual needs



shown with optional accessories

CDS-178

CENTRAL SUPPLY CART

OVERALL SIZE

66⅝"W x 26⅛"D x 69"H

FEATURES

- All welded stainless steel construction for one single piece frame
- Outboard plate-mounted 8" casters with diagonally mounted swivel locks provide stability
- Wrap around bumper mounted away from frame for easy cleaning
- Includes nine standard 5" deep tote boxes or up to four large 10" deep tote boxes
- Optional shelves and accessories available to configure the cart to suit individual needs





FEATURES

- All welded stainless steel construction
- All seams fully welded
- Frame tubing is 1¼", 18 gauge stainless steel
- Wrap around bumper mounted away from the frame for easy cleaning
- Shelves adjust every four inches
- Low maintenance
- Solid stainless steel pan bottom with reinforcing channel and sound control strip for added strength and quiet operation



CDS-147

OVERALL SIZE

39⅜"W x 26⅛"D x 42⅛"H

CDS-147-A

OVERALL SIZE

39⅜"W x 26⅛"D x 54⅛"H

- Ideal for case cart systems or general floor use
- Shelf size 35¼"W x 20⅜"D
- 5" swivel casters



Swivel Caster



CDS-148

OVERALL SIZE

50⅛"W x 26⅛"D x 58"H

- Suitable for general floor service, pharmacy or linen
- Shelf size 46⅜"W x 20⅜"D
- 8" plate casters with two diagonally mounted swivel locks



Plate Caster



CDS-149

OVERALL SIZE

60½"W x 26⅛"D x 66⅛"H

- Ideal cart for linen distribution, storage and floor nursing stations
- Shelf size 56⅜"W x 20⅜"D
- 8" plate casters with two diagonally mounted swivel locks



Plate Caster

ACCESSORIES



Two Way Tote Boxes

- 5" deep
- 10" deep



Adjustable Shelf Dividers

- 6" tall
- 10" tall



Wire Backs

- 2" x 3" openings
- 4" x 6" openings



Adjustable Retaining Rods

- Single style
- Double style



Stainless Steel Shelves

- Wire
- Solid





UTILITY CARTS



CDS-140

Light Weight Utility Cart

OVERALL SIZE

27½"W x 18½"D x 40"H

SHELF WEIGHT CAPACITY

Top 25 lbs.

Middle 25 lbs.

Bottom 50 lbs.

Maximum Load 100 lbs.

FEATURES

- Stainless-steel, Steril-Gard construction
- Lightweight 1" by 16-gauge tubular frame
- Three shelves welded to legs and turned up 1" at the edges
- Outboard stem mounted 4" casters for stability



CDS-140-B

Light Weight Utility Cart

OVERALL SIZE

32"W x 21¼"D x 40"H

SHELF WEIGHT CAPACITY

Top 25 lbs.

Middle 25 lbs.

Bottom 50 lbs.

Maximum Load 100 lbs.

FEATURES

- 5" rolling bumpers to protect walls and door jambs
- Stainless-steel, Steril-Gard construction
- Lightweight 1" by 16-gauge tubular frame
- Three shelves welded to legs and turned up 1" at the edges
- Outboard stem mounted 4" casters for stability



CDS-145

Utility Cart

OVERALL SIZE

41½"W x 22"D x 38¼"H

SHELF WEIGHT CAPACITY

Top 50 lbs.

Middle 50 lbs.

Bottom 100 lbs.

Maximum Load 200 lbs.

FEATURES

- "A-frame" design allows easier access from sides
- Stainless-steel, Steril-Gard construction
- Easy to clean
- Strong, fully welded construction
- Three sides turned up to retain items, front side turned down for easy access

optional Hand Grips available (CDS-145-HG)



NEW!

CDS-153-HD

Heavy-Duty Utility Cart

OVERALL SIZE

39¾"W x 27¾"D x 44¾"H

SHELF WEIGHT CAPACITY

Top 250 lbs.

Middle 250 lbs.

Bottom 400 lbs.

Maximum Load 900 lbs.

FEATURES

- Strong, fully welded
- Stainless-steel, Steril-Gard construction
- Easy to clean

optional brake casters available - contact customer service for details





MULTI-PURPOSE CARTS

MULTI-PURPOSE CARTS



shown with optional linen ring

CDS-153

CASE CART

OVERALL SIZE

36½"W x 26¾"D x 35½"H

FEATURES

- Optional linen rings to hold hamper or trash bags
- Linen ring locks in up position and swings down when not in use
- Equipped with rails on all three surfaces
- Adjustable stainless steel center shelf
- Wrap around bumper
- 5" diameter swivel plate casters



Swivel Plate Caster



shown with optional accessories

CDS-262

MULTI USE PROCEDURE CART

OVERALL SIZE

29⅝"W x 29⅝"D x 62"H

FEATURES

- Designed to be used as a surgical or obstetrical case cart.
- Multipurpose cart can also be used as a linen cart or as a clean/dirty utility cart
- May be ordered with optional shelves, wire backs, and other accessories
- Wrap around bumper
- 6" diameter swivel plate casters



Swivel Plate Caster



shown with optional accessories

CDS-270

UTILITY CART

OVERALL SIZE


50⅝"W x 26⅛"D x 65"H

FEATURES

- Designed to be used as a clean/dirty linen cart or as a large utility cart
- May be ordered with optional shelves, wire backs, and other accessories
- Wrap around bumper
- 6" diameter swivel plate casters with diagonally mounted swivel locks



Swivel Plate Caster

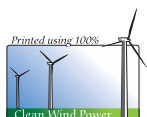






WIRE SHELVES

- Electro-polished stainless steel
- Easy to adjust or remove
- $\frac{3}{16}$ " diameter wire welded to $\frac{1}{4}$ " diameter wire support rods
- Maximum opening between rods: 1" x $6\frac{1}{2}$ "



SOLID SHELVES

- Steril-Gard construction
- Easy to adjust or remove
- 18-gauge stainless steel

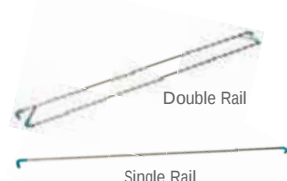

RETAINING RAILS

- Retaining rails to prevent accidental loss of materials
- Easily attached or removed
- Sound deadening rubber sleeve
- Polished stainless steel ($\frac{5}{16}$ " diameter)



SHELF DIVIDERS

- Adjustable
- Stainless steel construction for durability and ease of maintenance
- Available for solid or wire shelves in 6" or 10" height



5" TOTE BOX DRAWERS

- Constructed of impact-resistant, high-density polyethylene for durability
- Optional short or long tote box dividers available

DIMENSIONS:

17 $\frac{3}{4}$ "L x 14 $\frac{1}{2}$ "W x 5 $\frac{1}{8}$ "D

10" TOTE BOX DRAWERS

- Constructed of impact-resistant, high-density polyethylene for durability
- Optional short or long tote box dividers available

DIMENSIONS:

17 $\frac{3}{4}$ "L x 14 $\frac{1}{2}$ "W x 10"D

TOTE BOX HANGERS

- Hangs below both wire or solid shelving



Tote mounted with hanger

DISPOSABLE CART COVERS

- Heavy, clear plastic material
- Can serve as a trash bag after completion of the procedure
- Packed on convenient rolls
- Economical, single-use alternative to washable covers



HEAVY-DUTY COVERS

- Heavy duty, anti-bacterial STAPH-CHEK™ covers protect cart and its cargo
- Flame-resistant
- Clear vinyl front with heavy-duty nylon zippers or Velcro closures
- Washable, easy to maintain



WIRE BACKS

- Electro-polished stainless steel
- $\frac{3}{16}$ " diameter vertical wires, welded to $\frac{1}{4}$ " diameter horizontal rods
- Choice of 2" x 3" or 4" x 6" openings



AUTOMATIC COUPLERS

- Cart couplers are available for installation on all open carts
- Compatible with any major elevator or cart washer with automatic cart pickup system



STAINLESS STEEL CASTERS

- Ideal for carts that are cart washed daily
- All stainless steel material and bearings for a rust-free caster
- Two lubrication points for ease of maintenance



Stainless Steel Caster

SWIVEL LOCK

- Used for improved steering
- Recommended for long carts to assist in navigating long corridors.
- Disengage swivel lock for maneuverability in tight spaces
- Easy foot operation



Swivel Lock





CLOSED CART ACCESSORIES



WIRE SHELVES

- Electro-polished stainless steel
- Easy to adjust or remove
- $\frac{3}{16}$ " diameter wire welded to $\frac{1}{4}$ " diameter wire support rods
- Maximum opening between rods: 1" x 6 $\frac{1}{2}$ "
- Option of roll-out or stationary shelves



SOLID SHELVES

- Steril-Gard construction
- Easy to adjust or remove
- 18-gauge stainless steel
- Option of roll-out or stationary shelves



Positive Leaf Lock

DROP LEAF SHELVES

- Extends top of cart to increase work surface area
- Drain holes at end of leaf
- Features positive leaf lock



TOTE BOX DRAWERS

- Constructed of impact-resistant high density polyethylene for durability
- Available in 5" or 10" sizes
- Optional short or long tote box dividers available to section drawer into smaller compartments



TOTE BOX HANGERS

- Stainless steel hangers hold totes below top of the cart
- Attaches to either wire or solid shelving



CART IDENTIFICATION SYSTEMS

- For easy cart identification
- High-visibility Lexan material
- Cart washable
- Custom labels available upon request

CARD HOLDER

- Standard card sizes include:
 - 4" x 6"
 - 5" x 7"
 - 8 $\frac{1}{2}$ " x 11"
- Available in any size or orientation to suit individual needs



EXTERIOR SECURITY LATCH

- Prevents unintentional opening of cart doors
- Available with hole to fit customer provided tamper proof tags



Tamper Proof Tag

AUTOMATIC COUPLERS

- Cart couplers are available for installation on all closed carts
- Compatible with any major elevator or cart washer with automatic cart pickup system



STAINLESS STEEL CASTERS

- Ideal for carts that are cart washed daily
- All stainless steel material and bearings for a rust-free caster
- Two lubrication points for ease of maintenance



Stainless Steel Caster

SWIVEL LOCK

- Used for improved steering
- Recommended for long carts to assist in navigating long corridors.
- Disengage for maneuverability in tight spaces.
- Easy foot operation




Swivel Lock

BRAKE LOCK

- Prevents cart from rolling
- Easy foot operation



Brake Lock

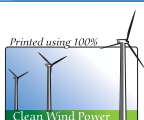



SPACE STATION



SPACE STATION

A specialty back table designed with additional room for large cases such as orthopedics, spinal fusions, neuro and craniotomies. The tuck-away cantilevered shelf holds multiple trays and is angled for clear observation of instruments.

The Space Station is easily rolled into position using 4" diameter, heavy-duty ball bearing brake/swivel casters. To lock in place, simply engage the single-pedal brake and swivel caster locks.

CDS-3060

Main Table
Top Shelf

Dimensions

30"D x 60"W x 34"H
20"D x 60"W

CDS-3072

Main Table
Top Shelf

Dimensions

30"D x 72"W x 34"H
20"D x 72"W



ONE PERSON OPERATION



Pneumatic assist enables one-person operation in raising or lowering the cantilevered shelf.

Pneumatic assist



Shelf is easily raised when extra space is required, or lowered into its storage position if not needed.

160 lbs.
Weight Capacity

230 lbs.
Weight Capacity

230 lbs.
Weight Capacity



Cantilevered shelf features a 10° slope and has an adjustable height from 49.9" to 52.9".

DISPOSABLE TABLE DRAPES


Single piece disposable sterile drapes are available to cover both surfaces of the sterile field.

CDS-3060-DD

12 Drapes per case

CDS-3072-DD

12 Drapes per case

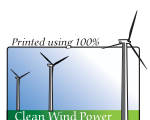

800.822.3501 • FAX 360.696.1700 • www.pedigo-usa.com



4000 SE Columbia Way • Vancouver, WA 98661



08/09 MM-60



INSTRUMENT
ASSEMBLY STATIONS

INSTRUMENT ASSEMBLY STATIONS



CDS-1000

MAIN TABLE TOP SURFACE

66"W x 30"D

HEIGHT ADJUSTABLE


36" to 47½"H

FEATURES

- Wrap-around work surface for easy access to all areas
- Electrically operated hydraulic height adjustment provides comfortable, ergonomic height for all staff members
- Three large pegboards for organizing instruments for storage and easy access
- Eight slide-out lower shelves for bins
- 3" diameter, heavy-duty locking casters for easy mobility
- Stainless steel construction
- Six removable stainless steel trays, recessed into the top work surface to hold small items
- Center drawer standard (keyboard tray option available)



Removable Stainless Steel Trays

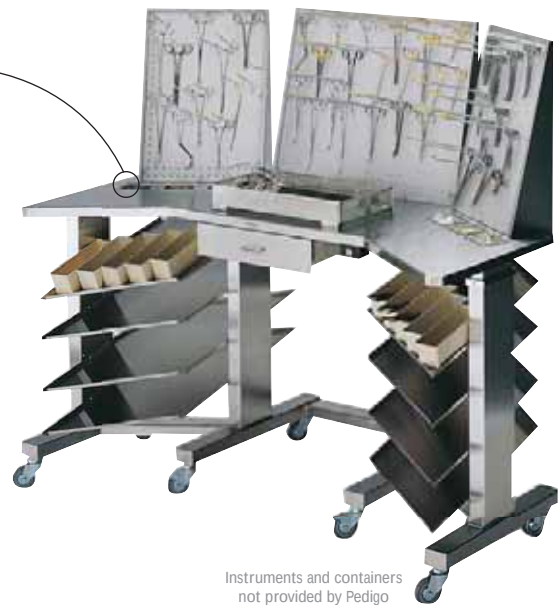

Instruments and containers not provided by Pedigo

CDS-1100

MAIN TABLE TOP SURFACE

72"W x 30"D

HEIGHT ADJUSTABLE

32" to 43¼"H

FEATURES

- Electrically operated hydraulic height adjustment provides comfortable, ergonomic height settings for all staff members.
- Shelves are positioned 18" and 30" above main work surface, and measure 72"W x 12"D with 5/8" front lip.
- Shelves also feature a 10° slope.
- Load capacity of 225 pounds on main work surface and 50 pounds on each upper shelf.
- 4" stainless steel brake/swivel lock casters with precision ball bearings for easy mobility.



OPTIONAL ACCESSORIES

1 Work Lamp CDS-1020

- Work Lamp arm has a 42" reach
- Features a 5" viewing lens
- 75% magnification

2 Outlet Strip CDS-1015

- Hospital grade
- Six outlets
- 15' foot cord

3 Keyboard Tray CDS-1032

- Measures 22¼"W x 10"D x 2¾"H
- Tray extends 10½" (retractable)
- Built-in mouse tray comes standard

4 Computer Enclosure CDS-1031

- Open back for air flow and easy access to computer cords.
- Enclosure can be customized to accommodate larger computer towers (max. load 30 lbs/13.60 kg)

5 Extension Arm CDS-1030

- Supports LCD monitors (max. load 25 lbs / 11.3 kg)
- For optimal viewing, extension arm features multi-directional adjustment
- Gas-strut and ball-joint system for smooth, stable operation

6 Storage Cart CDS-1010

- 4 drawer tote box cart for added storage fits under the workstation
- Totes are constructed of impact-resistant, high-density polyethylene and measure 17¾"L x 14½"W x 5"D





SPECIALTY OR TABLES

SPECIALTY OPERATING ROOM TABLES



SPACE STATION Specialty Back Table



A specialty back table designed with additional room for large cases such as orthopedics, spinal fusions, neuro and craniotomies. The tuck-away cantilevered shelf holds multiple trays and is angled for clear observation of instruments.

The Space Station is easily rolled into position using 4" diameter, heavy-duty ball bearing brake/swivel casters. To lock in place, simply engage the single-pedal brake and swivel caster locks.

- Cantilevered shelf for easy observation of instruments
- One-person operation in raising or lowering shelf

CDS-3060

Main Table
Top Shelf

Dimensions

30"D x 60"W x 34"H
20"D x 60"W

CDS-3072

Main Table
Top Shelf

Dimensions

30"D x 72"W x 34"H
20"D x 72"W

DISPOSABLE TABLE DRAPES

Single piece disposable sterile drapes are available to cover both surfaces of the sterile field.

CDS-3060-DD

12 Drapes per case

CDS-3072-DD

12 Drapes per case



OVER-OPERATING Specialty Table



Pedigo's Over-Operating Specialty Table is constructed of stainless steel and features easy, one-hand height adjustment.


- 3" total lock casters, with brakes and swivel locks
- Adjustable height from 39" to 54¼" high from floor
- 27³/₈" clearance between legs
- Optional removable extension leaf adds 10" to the length of the top
- Other sizes available upon request.

P-5182

Standard Top

Dimensions

34"W x 42"L

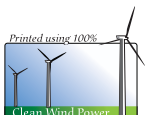



CENTRAL SUPPLY TABLES



shown with optional accessories

CDS-2472

OVERALL SIZE

72"W x 24"D x 34"H

FEATURES

- Heavy duty, all stainless steel construction
- Height adjustable floor glide
- For use in areas where a table larger than a standard utility or work table is needed

OPTIONAL ACCESSORIES

- Tote box drawers, 5" or 10" deep, or stainless steel drawers
- Full length lower shelf
- 4" swivel casters with brakes



shown with optional accessories

CDS-3260

OVERALL SIZE

60"W x 32"D x 34"H

FEATURES

- Heavy duty, all stainless steel construction
- Height adjustable floor guides

OPTIONAL ACCESSORIES

- Sterile Wrap frame, 56"W x 30½"H
- Tote box drawers, 5" or 10" deep, or stainless steel drawers
- Full length lower shelf
- 4" swivel casters with brakes



shown with optional accessories

CDS-3684

OVERALL SIZE

84"W x 36"D x 34"H

FEATURES

- Heavy duty, all stainless steel construction
- Height adjustable floor glide
- For use in areas where a table larger than a standard utility or work table is needed

OPTIONAL ACCESSORIES

- Tote box drawers, 5" or 10" deep, or stainless steel drawers
- Full length lower shelf
- 4" swivel casters with brakes



Caster with brake





INSTRUMENT TABLES



SG-86-SS

MID-SIZE TABLES WITH SHELF

DIMENSIONS

SG-82-SS	20"W x 16"D x 34"H
SG-84-SS	30"W x 16"D x 34"H
SG-86-SS	33"W x 18"D x 34"H
SG-88-SS	36"W x 20"D x 34"H



SG-87-SS

MID-SIZE TABLES WITHOUT SHELF

DIMENSIONS

SG-83-SS	20"W x 16"D x 34"H
SG-85-SS	30"W x 16"D x 34"H
SG-87-SS	33"W x 18"D x 34"H
SG-89-SS	36"W x 20"D x 34"H



SG-92-SS

LARGE TABLES WITH SHELF

DIMENSIONS

SG-90-SS	48"W x 20"D x 34"H
SG-92-SS	48"W x 24"D x 34"H
SG-94-SS	60"W x 24"D x 34"H
SG-96-SS	72"W x 24"D x 34"H
SG-98-SS	36"W x 24"D x 34"H



SG-93-SS

LARGE TABLES WITHOUT SHELF

DIMENSIONS

SG-91-SS	48"W x 20"D x 34"H
SG-93-SS	48"W x 24"D x 34"H
SG-95-SS	60"W x 24"D x 34"H
SG-97-SS	72"W x 24"D x 34"H

FEATURES

- Stainless steel Steril-Gard construction eliminates full sub-tops where undesirable substances can collect
- Lightweight but strong
- Attractive design
- Easy to clean
- Soundproofed reinforcing channels for quiet use
- 3" rubber, double ball bearing casters for silent mobility
- Available in many sizes

OPTIONS

- 3 or 4 sided top and shelf rails available for all tables
- Stainless steel drawers, or tote boxes may be mounted under top.
- Table height can be modified

Glides/ Rubber Feet/ Casters

- #15 Glides
- #13 Rubber Feet
- #06 Locking Casters

Tote Boxes

Totes are constructed of impact-resistant, high-density polyethylene for durability


- CDS-TB-5" Tote Box (17¾"L x 14½"W x 5⅞"D)
- CDS-TB-10" Tote Box (17¾"L x 14½"W x 10"D)

Tote Box Hangers

- SG-TBH-1 16"D SG Tables
- SG-TBH-2 18"D SG Tables
- SG-TBH-3 20"D SG Tables
- SG-TBH-4 24"D SG Tables

IV Pole for SG Tables

- P-7201-1 available on lower shelf models only.

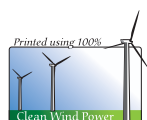



UTILITY TABLES



SG-80-SS

Utility Table with shelf, drawer and rail

OVERALL SIZE

20"W x 16"D x 34"H

DRAWER

15½"W x 12½"D x 3½"H



SG-80-A-SS

Utility Table with shelf, two drawers and rail

OVERALL SIZE

20"W x 16"D x 34"H

DRAWER

15½"W x 12½"D x 3½"H



SG-80-B-SS *no rail*

Utility Table with shelf and drawer

OVERALL SIZE

20"W x 16"D x 34"H

DRAWER

15½"W x 12½"D x 3½"H



SG-80-C-SS *no rail*

Utility Table with shelf and two drawers

OVERALL SIZE

20"W x 16"D x 34"H

DRAWER

15½"W x 12½"D x 3½"H



SG-81-SS

Utility Table with shelf and rail

OVERALL SIZE


20"W x 16"D x 34"H


FEATURES

- Convenient and easy to use
- Stainless steel, Steril-Gard construction eliminates full sub-tops
- Improves the control of undesirable substances
- Full extending drawer hangers (*except SG-81-SS*)
- Sound-proofed reinforcing channels for quiet use and added strength
- 3" double ball bearing, rubber casters for silent mobility
- 1½" high 3-sided top rail

OPTIONS

- Part No. 46-B Drawer lock with key
- Part No. 06 3" casters with brakes

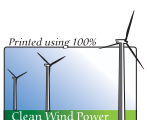

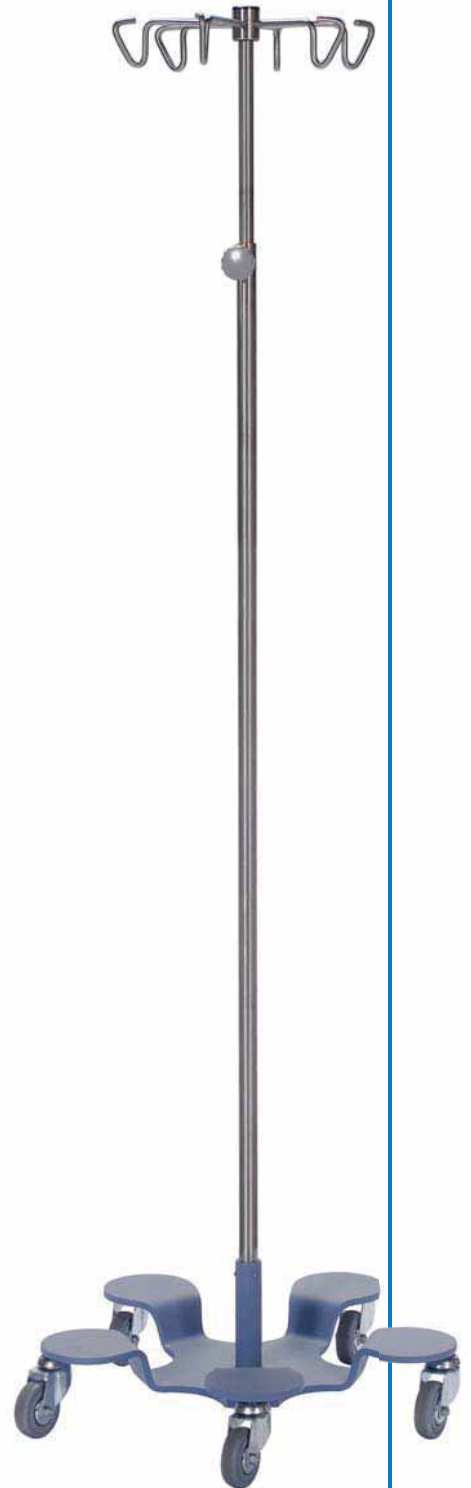

***i*STAND™**

INFUSION PUMP STAND



Pedigo prides itself on quality products manufactured here in the United States and built to last a lifetime.

The iSTAND™ will last a lifetime!

FEATURES

- Accommodates up to six infusion pumps
- Bag Drop Prevention Device
allows for controlled descent of pole
- Stainless steel pole and hooks
- Adjustable height from 69" to 124"
highest in the industry!
- Heavy duty 3" sealed casters
- Extremely durable powder coat finish base
- Pole sleeve welded securely to the base for reinforced stability

P-1080-6 iSTAND™



Pole

Stainless steel resists staining, corrosion and rusting.

Brass Bushing

Provides smooth operation while reducing galling and wear.

NEW! Infusion Pump Hoop

Pedigo's NEW Infusion Pump Hoop measures 29"H x 19"W and allows for multiple pumps to be installed. Made from one inch stainless steel tubing, the hoop allows for rotational and vertical adjustment.

P-1080-F Infusion Pump Hoop

Base

Low center of gravity design and 27 lbs. base offers improved stability.

Footprint - 23"

Small footprint reduces stumbles and falls while patient ambulates.

ACCESSORIES



Patient Handle

For easy mobility

P-3301 Stainless Steel
P-3302-BLK Black

includes Universal Clamp



Wire Basket

P-3201 Large, Stainless
P-3202 Large, Black
P-3211 Medium, Stainless
P-3212 Medium, Black

includes Universal Clamp



Infusion Pump Holder

Multiple pump holder

P-1080-C Infusion Pump Holder

includes Universal Clamp



Foley Hook

Multiple pump holder

P-1080-E Foley Hook

includes Universal Clamp



IV Set Organizer

Also holds secondary bags

P-1080-B Stainless Steel



Outlet Strip

Hospital grade

P-3401 Six outlet strip

includes Universal Clamp



O₂ Holder

Holds E-size oxygen tanks

P-3601 Stainless Steel

includes Universal Clamp




Universal Clamp

Required for mounting most iSTAND accessories

P-3500 Universal Clamp



Accessories above may be adapted to a majority of IV stands by utilizing the Universal Clamp.

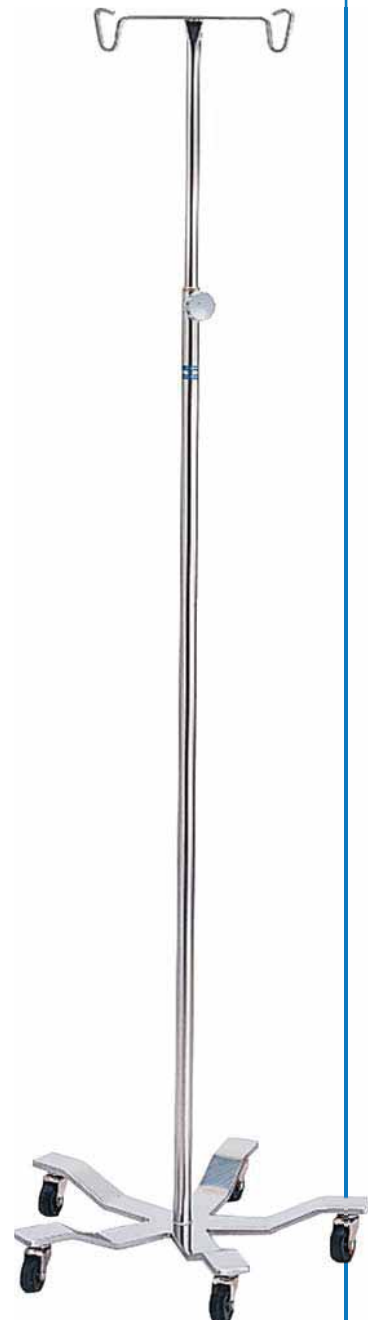

IV STANDS

IV STANDS



FEATURES

- Heavy, tip-resistant, all welded design
- 2" swivel casters for easy mobility
- 21" base on all stands
- 68" to 116" height adjustment
- Tube: O.D. = 1" I.D. = 3/4"
- Easy to clean surfaces

P - 1 0 7 2 - 2 *two hook*


P - 1 0 7 2 - 4 *four hook*

- Stainless steel
 - Hand operated, friction lock height Adjustment
 - Four leg base
 - 70 lbs. overall weight capacity
 - 10 lbs. per hook
(recommended 30 lbs. max. evenly distributed on top hooks for stability)
- ...see Safe Load Capacity

P - 7 6 *two hook*

P - 7 6 - 4 *four hook*

- Chrome
 - Hand operated, friction lock height Adjustment
 - Four leg base
 - 70 lbs. overall weight capacity
 - 10 lbs. per hook
(recommended 30 lbs. max. evenly distributed on top hooks for stability)
- ...see Safe Load Capacity

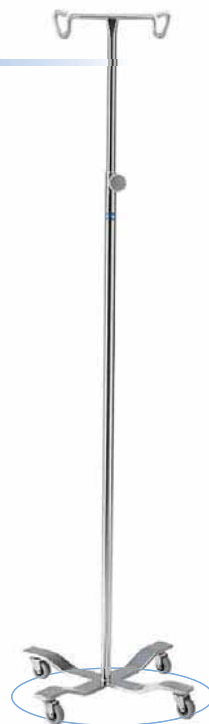

19" Caster Center Line



18 1/2" Caster Center Line

P - 1 5 7 2 - 2 *two hook*

P - 1 5 7 2 - 4 *four hook*

- Stainless steel
 - Hand operated, friction lock height adjustment
 - Five leg base
 - 70 lbs. overall weight capacity
 - 10 lbs. per hook
(recommended 30 lbs. max. evenly distributed on top hooks for stability)
- ...see Safe Load Capacity

P - 5 7 6 - 2 *two hook*

P - 5 7 6 - 2 *four hook*

- Chrome
 - Knock down design for shipping
 - Five leg base
 - 70 lbs. overall weight capacity
 - 10 lbs. per hook
(recommended 30 lbs. max. evenly distributed on top hooks for stability)
- ...see Safe Load Capacity

P - 1 0 7 6 - 2 *two hook*


P - 1 0 7 6 - 4 *four hook*

P - 1 0 7 6 - 6 *six hook*


- Stainless steel
 - Quick and easy one-hand foot pedal height adjustment
 - 70 lbs. overall weight capacity
 - 10 lbs. per hook
(recommended 30 lbs. max. evenly distributed on top hooks for stability)
- ...see Safe Load Capacity



Foot pedal height adjustment

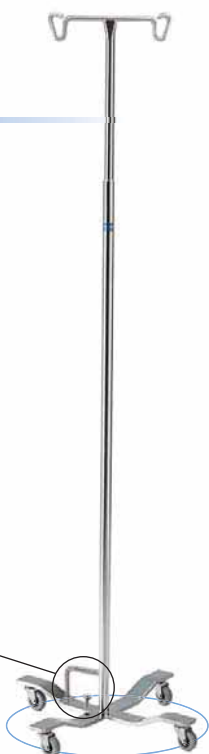

19" Caster Center Line

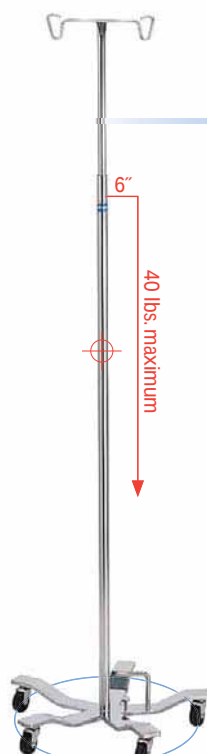

18 1/2" Caster Center Line

P - 1 5 7 6 - 2 *two hook*

P - 1 5 7 6 - 4 *four hook*

P - 1 5 7 6 - 6 *six hook*

- Stainless steel
 - Quick and easy one-hand foot pedal height adjustment
 - 70 lbs. overall weight capacity
 - 10 lbs. per hook
(recommended 30 lbs. max. evenly distributed on top hooks for stability)
- ...see Safe Load Capacity

Safe Load Capacity

For stability, a safe load on a Pedigo I.V. pole should not exceed 40 lbs. to one side and needs to be within 6" of base tube center line.





MAYO STANDS



P - 1 0 6 5 - S S Foot Operated

- Removable, extra large, 20" x 25" stainless steel tray for easy cleaning
- Convenient foot-pedal operated lock adjusts height from 36" to 62"
- Base designed to move under low clearance equipment with maximum stability
- All welded stainless steel construction
- 2" double ball bearing swivel casters



P - 1 0 6 6 - S S Mayo Stand

- Removable 16¼" x 21¼" stainless steel tray on stainless steel frame
- Convenient foot-pedal operated lock adjusts height from 36" to 62"
- Base designed to move under low clearance equipment with maximum stability
- All welded stainless steel construction
- 2" double ball bearing swivel casters



P - 1 0 6 9 - S S Foot Operated

- Removable, large, 16¼" x 21¼" stainless steel tray for easy cleaning
- Convenient foot-pedal operated lock adjusts height from 36" to 62"
- Three leg base design for increased mobility
- All welded stainless steel construction
- 2" double ball bearing swivel casters

P - 1 0 6 8 - S S Foot Operated

Tray Dimensions 12⅝" x 19⅛"

P - 1 0 6 8 - A - S S Hand Operated

Tray Dimensions 12⅝" x 19⅛"

P - 1 0 6 6 - A - S S Hand Operated

Tray Dimensions 16¼" x 21¼"

- Adjust height from 32" to 53" - friction lock

P - 1 0 6 9 - A - S S Hand Operated

Tray Dimensions 16¼" x 21¼"



P - 6 5 Medical Tray Unit

- Removable 19-1/8" L x 12-5/8" W stainless steel tray on chrome plated frame
- Hand operated, friction lock adjusts height from 30" to 40"
- 2" double ball bearing swivel casters

P - 6 6 Mayo Stand

- Removable 19-1/8" L x 12-5/8" W stainless steel tray on powder-coated frame
- Hand operated, friction lock adjusts height from 32" to 53"
- 2" double ball bearing swivel casters





BASIN STANDS &
KICK BUCKETS

BASIN STANDS



P-78 chrome basin stand

single basin stand

P - 7 8

Chrome

14"D x 14"W x 34"H

- One removable 13 $\frac{1}{4}$ " diameter, 7 quart stainless steel basin
- Stainless steel construction

P - 1 0 7 8 - S S

Stainless steel

14"D x 14"W x 34"H

FEATURES

- Elevated ring basin bracket
- Open frame design allows access from all sides
- 2" diameter double ball bearing swivel casters
- Mobile and easy to clean



P-78 chrome basin stand

single basin stand with shelf

P - 1 0 7 8 - W / S - S S

Stainless steel

14"D x 14"W x 34"H

- One removable 13 $\frac{1}{4}$ " diameter, 7 quart stainless steel basin
- Stainless steel construction

double basin stand with shelf

P - 1 0 7 9 - W / S - S S

Stainless steel

14"D x 28"W x 34"H



P-79 chrome basin stand

double basin stand

P - 7 9

Chrome

14"D x 28"W x 34"H

- Two removable 13 $\frac{1}{4}$ " diameter, 7 quart stainless steel basin
- Stainless steel construction

double basin stand

P - 1 0 7 9 - S S

Stainless steel

14"D x 28"W x 34"H

KICK BUCKETS



kick bucket

P - 1 0 2 0 - S S

Stainless steel

Overall size
13 $\frac{1}{2}$ " diameter x 14"H

- Removable 12 quart, stainless steel basin

FEATURES

- Balanced, tip-resistant
- Electro-polished, stainless steel frame
- Perimeter bumper to protect walls
- 2" swivel casters for easy mobility




sponge receptacle

P - 1 0 2 5 - S S

Stainless steel

Overall size
6 $\frac{1}{4}$ " diameter x 8"H

- Removable 8 $\frac{1}{2}$ quart, stainless steel basin

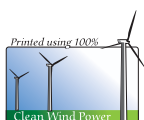

800.822.3501 • FAX 360.696.1700 • www.pedigo-usa.com



4000 SE Columbia Way • Vancouver, WA 98661



08/09 MM-69



LINEN HAMPERS



L I N E N H A M P E R S

Pedigo's newly designed linen hampers come standard with Quiet-Close™ lid dampener for improved safety and added control of airborne bacteria and odors.

Linen hampers are constructed of heavy gauge steel and feature a triangular shape for efficient use of space. Double ball bearing 3" swivel casters make for easy mobility. Strong, heavy-gauge shelf keeps bag from touching floor and simplifies bag removal.



P-120-L and P-121-L



P-1120-L-SS and P-1121-L-SS



P-120 and P-121



P-1120-SS and P-1121-SS

P - 1 2 0 - L

with Quiet-Close™ Lid

Dimensions

32"H with a 57¼" perimeter

- Powder coated steel
- Foot operated
- Accommodates 18" bags

P - 1 1 2 0 - L - S S

with Quiet-Close™ Lid

Dimensions

32"H with a 57¼" perimeter

- Stainless steel
- Includes rolling bumpers
- Foot operated
- Accommodates 18" bags

P - 1 2 0

open linen hamper

Dimensions

32"H with a 57¼" perimeter

- Powder coated steel
- Accommodates 18" bags

P - 1 1 2 0 - S S

open linen hamper

Dimensions

32"H with a 57¼" perimeter

- Stainless steel
- Rolling bumpers
- Accommodates 18" bags

P - 1 2 1 - L

with Quiet-Close™ Lid

Dimensions

35"H with a 78¼" perimeter

- Powder coated steel
- Foot operated
- Accommodates 25" bags

P - 1 1 2 1 - L - S S

with Quiet-Close™ Lid

Dimensions

35"H with a 78⁵/₈" perimeter

- Stainless steel
- Includes rolling bumpers
- Foot operated
- Accommodates 25" bags

P - 1 2 1

open linen hamper

Dimensions

35"H with a 78⁵/₈" perimeter

- Powder coated steel
- Accommodates 25" bags

P - 1 1 2 1 - S S

open linen hamper

Dimensions

35"H with a 78⁵/₈" perimeter

- Stainless steel
- Includes rolling bumpers
- Accommodates 25" bags

Quiet-Close™ Lid Damper

Pedigo's Quiet-Close™ hydraulic lid dampener slows the closure of the lid reducing the spread of bacteria and protects fingers and hands from being pinched.

- Foot operated
- Eliminates slamming
- Quiet operation
- Tested for 50,000 cycles



25" Bag Hamper Opening



18" Bag Hamper Opening



* approximately the size of a basketball





BASSINETS





shown with optional
bassinet and pad

P - 1110 - S S BASSINET STAND WITH DRAWER

FRAME SIZE

18¹/₈"W x 31"L x 37³/₄"H

Height with bassinet basket - 41"H

P - 1110

- Same as P-1110-SS except frame is powder coated steel

FEATURES

- Design allows close approach from both sides
- Basket can be raised at either end
- 3" double ball bearing swivel casters
- Easy to clean
- Stainless steel construction



shown with optional
bassinet and pad

P - 1110 - A - S S BASSINET STAND WITH DRAWER

FRAME SIZE

18¹/₈"W x 31"D x 37³/₄"H

Height with bassinet basket: 41"H

P - 1110 - A

- Same as P-1110-A-SS except frame is powder coated steel

FEATURES

- End mounted drawer measures: 17"W x 15¹/₂"D x 6¹/₂"H
- Design allows close approach from both sides
- Basket can be raised at either end
- 3" double ball bearing swivel casters
- Easy to clean
- Stainless steel construction



shown with optional
bassinet and pad

P - 1110 - B - S S BASSINET STAND WITH DRAWER

FRAME SIZE

18¹/₈"W x 31"D x 37³/₄"H

- Height with bassinet basket: 41"H

FEATURES

- Side mounted drawer measures: 21"W x 15¹/₄"D x 3³/₈"H
- Design allows close approach from both sides
- Basket can be raised at either end
- 3" double ball bearing swivel casters
- Easy to clean
- Stainless steel construction



shown with optional
bassinet and pad

P - 1114 - S S BASSINET STAND WITH DRAWERS

DIMENSIONS

32"W x 20"D x 33"H

Height with bassinet basket: 41"H

FEATURES

- Design consistent with rooming in techniques
- Two 10" deep tote boxes accessible from either side
- 4" double ball bearing swivel casters
- Easy to clean
- Stainless steel construction



shown with optional
bassinet and pad

P - 1115 - A - S S BASSINET CABINET

DIMENSIONS

30"L x 17¹/₂"W x 35"H

Drawer frame with removable pan

30³/₄"L x 16"W x 9¹/₂"H

P - 1115 - S S

Same as P-1115-A-SS except without doors

FEATURES

- Bassinet carriage slides back to allow working area
- Friction lock secures the bassinet in both open and closed positions
- Removable internal shelf, positioned in middle left compartment
- 4" double ball bearing swivel casters
- Easy to clean
- Stainless steel construction

OPTIONS


BASSINET BASKET (#53) BASSINET PAD (#54)

Basket measures

16"W x 30³/₄"D x 9⁵/₈"H

Pad measures

12¹/₂"W x 26¹/₂"D x 1" Thick

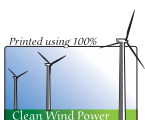



LABOR & DELIVERY

featuring
**RADIANT HEAT
TECHNOLOGY**



P - 2 0 1 0

BLANKET WARMING CABINET

COMPARTMENT SIZE

- 2.3 cubic feet
- 13¼"W x 21½"D x 14"H

OVERALL SIZE

16½"W x 25⅞"D x 20⅝"H

FEATURES

- No fans or other moving parts
- No lint traps
- 4 digit LED display
- Temperature adjusts to 200°F
- Optional casters or legs
- Safety coated interior
- 125 VAC



P - 1 1 1 0 - A - S S

BASSINET STAND WITH DRAWER

FRAME SIZE

18½"W x 31"D x 37¾"H

HEIGHT WITH BASSINET BASKET

41"H

OPTIONAL BASSINET BASKET (#53)

16"W x 30¼"D x 9⅝"H

OPTIONAL PAD (#54)

12½"W x 26½"D x 1" Thick

FEATURES

- Design allows close approach from both sides
- Basket can be raised at either end
- 3" double ball bearing swivel casters
- Drawer measures 17"W x 15½"D x 6½"H

Also available with side mounted drawer

Model P-1110-B-SS:

Drawer 21"W x 15¼"D x 3⅜"H

S T - 1 8 3 3 - S S

PEDIATRIC SCALE TABLE

OVERALL SIZE

33"W x 18"D x 35½"H

FEATURES

- Designed for use with common electronic scales
- Oversize 4" casters for easy movement
- One brake caster for stability

C D S - 2 5 6 - L D

SURGICAL CASE CART

CABINET SIZE

23"W x 24¼"D x 21"H

OVERALL SIZE

28⅝"W x 28⅝"D x 30"H

ACCESSORIES

- Wire shelves
- Solid shelves
- Roll-out shelves
- Drop-leaf shelves
- Card holders
- Stainless steel casters
- 5" and 10" tote boxes



shown with optional drop leaf shelf



ERGO

Task Chair



SEATING



P-521-GS

- 14" cushioned seat with back
- 23½" diameter aluminum base
- Pneumatic height adjustment
- 2" dual wheel casters
- 18¾" to 23¾" height adjustment



P-526-GS

- 14" round cushioned seat
- 23½" diameter aluminum base
- Pneumatic height adjustment
- 2" dual wheel casters
- 19¼" to 23¾" height adjustment



P-536

- 14" round cushioned seat
- 21½" diameter chrome base
- Screw shaft height adjustment
- Easy rolling 2½" casters
- 17¾" to 26⅛" height adjustment



P-536-GS

- 14" cushioned seat
- 21½" diameter chrome base
- Pneumatic height adjustment
- Easy rolling 2½" casters
- 18" to 23¼" height adjustment



P-551

- 16" x 17" contoured seat
- Adjustable back support
- Screw shaft height adjustment
- Easy rolling 2½" casters
- 20⅝" to 29" height adjustment



P-551-GS

- 16" x 17" contoured seat
- Adjustable back support
- Pneumatic height adjustment
- Easy rolling 2½" casters
- 18½" to 24" height adjustment



P-555-GS

- 16" x 17" contoured seat
- Adjustable back support
- Pneumatic height adjustment
- Easy rolling 2½" casters
- 25⅝" to 33½" height adjustment

*All Pedigo seating
has a weight capacity
rating of 300 lbs.*

STAINLESS STEEL



P-1037-SS

- 14" round stainless steel contoured seat
- Height adjustment from 20 $\frac{1}{8}$ " to 32 $\frac{3}{8}$ "
- Rubber floor tips



P-1038-SS

- 14" round cushioned seat
- Height adjustment from 22" to 34 $\frac{1}{2}$ "
- Rubber floor tips



P-1039-SS

- Upholstered spring back
- 14" round cushioned seat
- Height adjustment from 22 $\frac{3}{4}$ " to 34 $\frac{1}{2}$ "
- Rubber floor tips



P-1040-SS

- Upholstered spring back
- Cushioned square seat measures 14" x 15"
- Height adjustment from 26" to 37 $\frac{1}{2}$ "
- 3" double ball bearing swivel casters

P-1037-WIC-SS

- 2" double ball bearing swivel casters
- Height adjustment from 22 $\frac{1}{2}$ " to 34 $\frac{5}{8}$ "

P-1038-WIC-SS

- 2" double ball bearing swivel casters
- Height adjustment from 24 $\frac{1}{4}$ " to 36 $\frac{3}{4}$ "

P-1039-WIC-SS

- 2" double ball bearing swivel casters
- 25" to 36 $\frac{3}{4}$ " height adjustment

FOOT-OPERATED

P-527-GS

- Pneumatic height adjustment
- Adjustable back support
- 14 $\frac{1}{2}$ " round fully padded seat
- 19 $\frac{1}{4}$ " to 24 $\frac{1}{2}$ " height adjustment
- 23 $\frac{1}{2}$ " diameter aluminum base
- Five-casters




P-528-GS

- Pneumatic height adjustment
- 14 $\frac{1}{2}$ " round fully padded seat
- 19 $\frac{1}{4}$ " to 24 $\frac{1}{2}$ " height adjustment
- 23 $\frac{1}{2}$ " diameter aluminum base
- Five-casters



PEDIGO[®] Royal

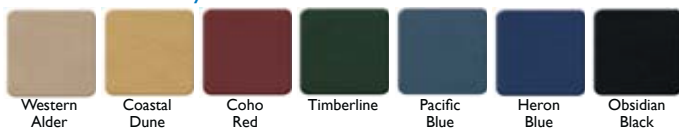

Environment friendly, Pedigo Royal upholstery is Latex-free, PVC-free, California TB-117 and TB-133 approved.



✓ LATEX-FREE ✓ PVC-FREE



PEDIGO[®] Classic



Cost effective, durable, vinyl based fabrics. Pedigo Classic upholstery is California TB-117 approved.



✓ LATEX-FREE



Above color samples are representations of upholstery materials. Actual colors may vary slightly. Material samples available upon request. You may also view color samples on our website at www.pedigo-usa.com.

P-35

- 14" cushioned seat
- Easy rolling 2" casters
- Friction lock height adjustment of 18" to 28 $\frac{1}{8}$ "

P-36

- Screw shaft height adjustment of 17 $\frac{7}{8}$ " to 26 $\frac{1}{4}$ "

P-38

- 14" cushioned seat
- Screw shaft height adjustment of 22" to 34 $\frac{1}{2}$ "
- Rubber floor tips

P-38-WIC

- 2" double ball bearing swivel casters
- Screw shaft height adjustment of 24 $\frac{1}{2}$ " to 36 $\frac{3}{4}$ "

P-39

- Upholstered spring back
- 14" round cushioned seat
- Screw shaft height adjustment of 22 $\frac{3}{4}$ " to 34 $\frac{1}{2}$ "
- Rubber floor tips

P-39-WIC

- 2" double ball bearing swivel casters
- Screw shaft height adjustment of 25" to 36 $\frac{3}{4}$ "

P-50

- Friction lock height adjustment of 18 $\frac{1}{2}$ " to 28 $\frac{1}{8}$ "
- Upholstered spring back
- 14" x 15" square cushioned seat

P-51

- Screw shaft height adjustment of 18 $\frac{5}{8}$ " to 26 $\frac{1}{4}$ "

P-50-R

- Friction lock height adjustment of 18 $\frac{1}{2}$ " to 28 $\frac{1}{8}$ "
- Upholstered spring back
- 14" round cushioned seat

P-51-R

- Screw shaft height adjustment of 18 $\frac{5}{8}$ " to 26 $\frac{1}{4}$ "

P-55 LAB STOOL

- Upholstered spring back
- 14" x 15" square cushioned seat
- Screw shaft height adjustment of 21 $\frac{5}{8}$ " to 33 $\frac{5}{8}$ "



ERGO TASK CHAIR

PEDIGO PROUDLY INTRODUCES

Our **NEW** ergonomically designed task chairs, the T-580, T-581 and T-582. All chairs feature pneumatic height adjustment, dual wheel casters, 5-star base, 300 lbs. weight capacity rating, integrated lumbar support, and independent backrest/seat pan height tilt controls.

✓ **LATEX-FREE**

✓ **PVC-FREE**



All task chair upholstery and plastic components meet California TB-117 and TB-133 flammability standards.



All task chair components meet the testing criteria for the American National Standard for Office Furnishings (ANSI/BIFMA X5.1-2002).

NEW!



T-580

- Adjustable arm rests
- 18" x 18½" contoured cushioned seat
- Pneumatic height adjustment of 17½" to 23¼"
- Easy rolling 2" casters
- 25" diameter footprint

T-581

- 18" x 18½" contoured cushioned seat
- Pneumatic height adjustment of 17½" to 23¼"
- Easy rolling 2" casters
- 25" diameter footprint

T-582

- 18" x 18½" contoured cushioned seat
- Pneumatic height adjustment of 17½" to 23¼"
- Easy rolling 2" casters
- 23" diameter footprint

HYDRAULIC SURGEON CHAIR

Hydraulic Surgeon Chairs are ideal for surgery, anesthesiology, ENT and reconstructive procedures. Engineered and constructed using our STERIL-GARD philosophy. STERIL-GARD is the concept and process of manufacturing products to minimize or eliminate any areas where dirt, fluid and bacterial debris can become entrapped. Therefore all Pedigo products are easier to clean, longer lasting, and most importantly - *reduce the transmission of bacteria*.

P-7000

- Foot-operated hydraulic height adjustment
- Single pedal height adjustment of 20 $\frac{3}{4}$ " to 28 $\frac{1}{4}$ "
- Specially contoured seat for added comfort
- Sensafoam body contouring seat and back
- Seat features tilt adjustment and backrest
- Height adjustable and articulating arm/wrist rests
- 21" diameter base with high grade 3" casters
- Chrome-plated 6-caster base
- Three locking casters



Locking Casters



P-6000


- Foot-operated hydraulic height adjustment
- Single pedal height adjustment from 22 $\frac{3}{4}$ " to 30 $\frac{1}{4}$ "
- Cushioned 16" diameter seat
- Seat features tilt adjustment and backrest
- 21" diameter base
- Powder-coated 5 caster base
- High grade 3" casters
- Three locking casters

P-6000 shown
without back rest



P-6000 OPTIONS

- Procedure Rest (Part No. P-6000-PR)
- Arm Rests (Part No. P-6000-D)

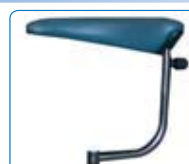

Procedure Rest



Arm Rests

P-6000 with Procedure Rest

- Rotates a full 360 degrees
- Upholstered for comfort
- Special density foam provides comfort and support during long procedures
- Cushioned 16" diameter seat
- Height adjustment from 24 $\frac{3}{8}$ " to 29 $\frac{1}{2}$ "



P-526 with Procedure Rest

- Rotates a full 360 degrees
- Upholstered for comfort
- Special density foam provides comfort and support during long procedures
- 14 $\frac{1}{2}$ " round fully padded seat
- Height adjustment from 16 $\frac{3}{4}$ " to 23 $\frac{3}{4}$ "

Part No. P-526-PR and P-526-GS

Part No. P-6000-PR and P-6000

* Seat tilt not available with procedure rest





PHYSICIAN
OFFICE EQUIPMENT



P-10 Footstool
300 lbs. weight capacity

P-10

Footstool • Chrome

- 16"W x 12"D x 8 1/4"H
- Strong, all-welded steel
- Rubber top mat and floor tips
- Optional handrail



P-66

Mayo Stand

- 13" x 19" stainless steel tray
- 32" x 53" height adjustment
- Easy 2" rolling casters
- lightweight yet strong & durable



SG-80-SS

Utility Table

OVERALL SIZE

20"W x 16"D x 34"H

- All stainless steel construction
- With guard rail and lower shelf

DRAWER

15 1/2"W x 12 1/2"D x 3 1/2"H



P-521-GS

Exam Stool

- 14" cushioned seat with back
- 24" diameter aluminum base
- Pneumatic height adjustment
- 2" dual wheel casters
- 19" to 25 1/2" height adjustment

Royal Upholstery



Environment friendly, Pedigo Royal upholstery is California TB-117 and TB-133 approved.



Classic Upholstery



Cost effective, durable, vinyl based fabrics. Pedigo Classic upholstery is California TB-117 approved.



manufacturing lifetime value

- Commitment to superior customer service and satisfaction
- Designed for ease of cleaning
- Large inventory for fast delivery
- Custom upholstery colors available to match every decor
- Pedigo items are found in physician offices and medical office buildings across the country
- Handcrafted & manufactured in the USA



other related items

- IV stands
- Basin stands
- Hampers
- Ergo Task Chair

Numerous stool models available



P-36

Exam Stool

- 14" cushioned seat
- Easy rolling 2" casters
- Screw shaft adjustment
- Height adjustment 17" to 27"

Above color samples are representations of upholstery materials. Actual colors may vary slightly. Material samples available upon request. All Pedigo upholstery features the protective PreFixx® coating for ease of cleaning and proven resistance to stains and abrasion. For more information and care instructions, please visit our website at www.pedigo-usa.com.





FEATURES

- Strong, all-welded steel
- Legs constructed of 1" O.D. 16-gauge tubular steel
- Handrail constructed of 1" O.D. tubular steel
- Tip resistant, rugged, and easy to maintain
- Corrugated, black, rubber top surface
- Rubber floor tips
- All carts built to Steril-Gard specifications



P-10 Footstool
300 lbs. weight capacity

footstool

P - 1 0

Chrome
12"D x 16"W x 8¼"H

P - 1 1

Chrome
12"D x 30"W x 8"H

P - 1 0 1 0 - S S

Stainless steel
12"D x 16"W x 8¼"H

P - 1 0 1 1 - S S

Stainless steel
12"D x 30"W x 5"H

P - 1 0 1 2 - S S

Stainless steel
12"D x 30"W x 8"H

with handrail

P - 1 0 - A

Chrome
12"D x 16"W x 8¼"H
• 38½"H with handrail

P - 1 0 1 0 - A - S S

Stainless steel
12"D x 16"W x 8¼"H
• 38½"H with handrail



P-10-A Footstool
300 lbs. weight capacity



P-18 Footstool
300 lbs. weight capacity

double step stool

P - 1 8

Chrome
12"D x 16"W x 16"H
• First step 8"H
• Second step 16"H

double step with handrail

P - 1 8 - A

Chrome
12"D x 16"W x 47"H with handrail
• First step 8"H
• Second step 16"H

stackable footstool

P - 1 0 1 5

Stainless steel
14"D x 19"W x 6"H


- Safely stack for increased height and interlocked to create larger platforms
- Slip resistant top
- Heavy gauge legs and handles

transport cart

P - 1 0 1 5 - T C

Stainless steel
12"D x 16"W x 8"H

- Transport cart holds up to six P-1015 for storage or mobility
- 3" diameter swivel casters

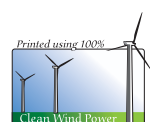

Pedigo proudly introduces our **NEW ergonomically designed anesthesia chair** the T-583. This chair features pneumatic height adjustment, dual wheel casters, 5-star base, integrated lumbar support and independent backrest/seat pan height tilt controls.


T-583 Ergo Task Anesthesia Chair



Environment friendly, Pedigo Royal upholstery is Latex-free, PVC-free, California TB-117 and TB-133 approved.



Seating Controls

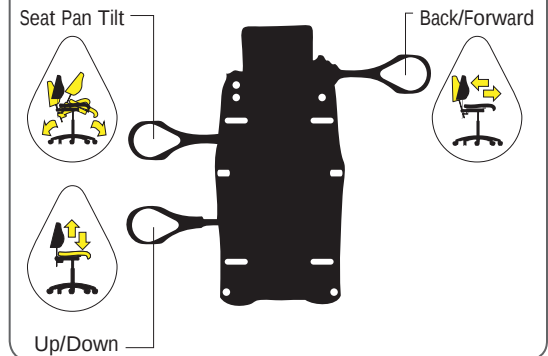



ivNow®

individual warming pods

Modular Configuration

Wall mount: 1-6 bag capacity

Countertop: 1-3 bag capacity

Pole mount: 1-3 bag capacity

Accommodates

0.5l, 1l, 2l and 3l bags

Anodized aluminum heat plate

Lightweight, medical grade plastic housing

Intelligent bag sensors

Dual temperature sensors

Heats to 104°F (40°C) in 30 minutes or less
(slightly longer heating time for 3l bags)



Fluid Warming Pods

Patent WO 2010/030746
Canada Patent 2692640



Shown with Pedestal Stand



Shown with Tilt Stand



ivNow-1

The ivNow-1 may be placed directly on the countertop or fitted with an optional countertop pedestal stand for additional stability.

P-2301-A

P-2301-B

Mount

- Countertop
- Equipment pole stand
- Wall using brackets

Capacity

One (1) bag, up to 3l *

ivNow-2

The ivNow-2 may be placed directly on the countertop or fitted with an optional countertop pedestal stand for additional stability.

P-2302-A

P-2302-B

Mount

- Countertop
- Equipment pole stand
- Wall using brackets

Capacity

Two (2) bags, up to 3l *

ivNow-3

The ivNow-3 unit requires either a mobile stand mount, wall mount bracket, tilt stand, or countertop pedestal stand prior to use.

P-2303-A

P-2303-B

Mount

- Countertop
- Equipment pole stand
- Wall using brackets

Capacity

Three (3) bags, up to 3l *

ivNow-4 • ivNow-5 • ivNow-6

ivNow wall mounted units require either the horizontal or vertical bracket to attached to a wall. See chart below for variable configurations.

P-2304-A

P-2305-A

P-2306-A

Mount

- Wall using brackets

Capacity

ivNow-4: Four (4) bags, up to 3l

ivNow-5: Five (5) bags, up to 3l

ivNow-6: Six (6) bags, up to 3l



P-1086-6 iSTAND™ shown above with accessories and P-2301-A ivNow

Accessories and Mounting Options

ivNow Units		Capacity (liters)	Number of Bags	Cord - Bottom (wall or mobile)	Cord - Back (countertop)	Pedestal Mount	
ivNow-1	P-2301-A	0.5l - 3l	1	•			
	P-2301-B	0.5l - 3l	1		•		P-2300-optic
ivNow-2	P-2302-A	0.5l - 3l	2	•			
	P-2302-B	0.5l - 3l	2		•		P-2300-optic
ivNow-3	P-2303-A	0.5l - 3l	3	•			
	P-2303-B	0.5l - 3l	3		•	P-2300-PM required	
ivNow-4	P-2304-A	0.5l - 3l	4	•			N/A
ivNow-5	P-2305-A	0.5l - 3l	5	•			N/A
ivNow-6	P-2306-A	0.5l - 3l	6	•			N/A


OPTIONS

Modular Thermal Fluid Warming

ivNow modular fluid warmers put warmed fluid at your fingertips where you need it, when you need it. ivNow fluid warmers are small, convenient units that can be placed in a number of different environments and combinations. Countertop, wall mounted and mobile units are all available depending on your particular requirements.

ivNow quickly warms and maintains safe temperatures of intravenous and irrigation fluids. The specially contoured warming surface cradles a single solution bag (0.5l - 3l). Each cavity is controlled by an individual control and L.E.D. display. A sensor in the heating plate detects the presence of a bag and engages the heating mechanism to quickly warm the fluid. Two temperature sensors continuously read the temperature of the fluid bag to engage the heater as necessary to maintain the fluid temperature within $\pm 0.2^{\circ}\text{C}$ ($\pm 0.3^{\circ}\text{F}$) of the set point temperature.

ivNow warmers ensure your fluids are kept at the proper temperature and will always be immediately available at or close to body temperature rather than room temperature without the added expense and set-up time required with in-line warming systems.

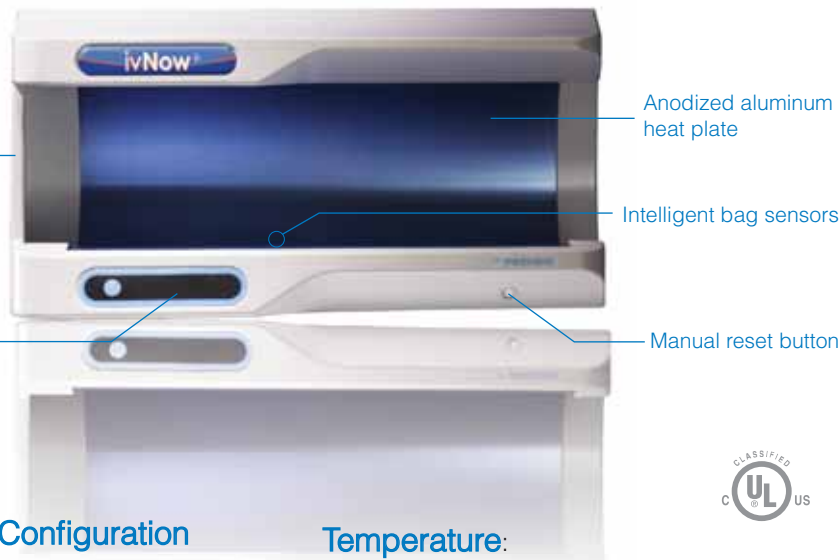

Modular Configuration

Wall mount:

One to six (1-6) bag capacity

Countertop:

One to three (1-3) bag capacity

Pole mount:

One to three (1-3) bag capacity

Warranty

1 year parts and labor

Temperature:

- Heats to 104°F (40°C) in 30 minutes or less (slightly longer heating time is required for 3L bags)
- Accuracy of $\pm 0.2^{\circ}\text{C}$ ($\pm 0.3^{\circ}\text{F}$)

Electronic Control

Monitors the length of time the bag has been held at temperature, displayed additionally by a status button. The control alerts users when a fluid bag has been held at temperature longer than 14 days.

	Mobile Stand Mount	Tilt Kit (required for 3l bags)	Wall Mount Kit (vertical)	Wall Mount Kit (horizontal)	Wall Mount Tilt Kit (horizontal)
	P-2300-MK		P-2300-VM	P-2300-HWM1	P-2300-WMT1†
0-PM onal		P-2300-CT1*			
	P-2300-MK		P-2300-VM	P-2300-HWM1	P-2300-WMT2†
0-PM onal		P-2300-CT2*			
	P-2300-MK		P-2300-VM	P-2300-HWM3	P-2300-WMT3†
		P-2300-CT3*			
A	P-2300-MK		P-2300-VM	P-2300-HWM3	N/A
A	P-2300-MK		P-2300-VM	P-2300-HWM3	N/A
A	P-2300-MK		P-2300-VM	P-2300-HWM3	N/A

† P-2300-WMT2

Wall Mount Tilt Kit requires Horizontal Wall Mount Kit (P/N P-2300-HWM1).

† P-2300-WMT3

Wall Mount Tilt Kit requires Horizontal Wall Mount Kit (P/N P-2300-HWM3).

* Tilt Kits may be configured using the Pedestal Mount Stand though not required.

NOTE:

Tilt kit is not recommended for 0.5l bags. Small bags may not activate the temperature power sensor when placed in the pod and tilted back.

Why buy ivNow?



It's easy to use

- Simply place an IV bag in the ivNow to activate the unit.
- Fluid is automatically warmed to 104°F (40°C) in less than 30 minutes.

It's good to the environment

- Eliminates the need for disposables associated with in-line fluid warmers, dramatically reducing hazardous disposable waste.
- Activates only when a fluid bag is on the warming tray and turns off when the bag is removed.
- Reduces electrical costs.

It saves money

- Save thousands of dollars annually on captive disposables that in-line warmers require.
- Eliminate time spent setting up and handling in-line intravenous warming disposables.
- With these savings, the ivNow can pay for itself in months.

It fits anywhere

- Innovative design allows ivNow to be kept conveniently near the patient on a countertop, wall, or pole mount.
- Soft colors and rounded corners blend with any modern healthcare décor



Scanning the QR code to the right with your mobile device will allow for playback of a video demonstration & features of the ivNow warming pods.





Elite Series

Blanket Warming Cabinets

Temperature Range

98°F - 180°F (37°C - 82°C)

Energy Consumption

Efficient, balanced control allows for low energy use and minimal heat loss.

Radiant Heat Technology

Utilizes multiple zone heating that monitors surfaces for safe temperature ranges.

Interior Lighting

Casts a comforting blue glow with two (2) different intensity settings.

Epoxy-Coated

interior surfaces

Three cabinet sizes available.

featuring

**RADIANT HEAT
TECHNOLOGY**



P-2220 • P-2210 • P-2205

Specifications:

- Single-chamber warming cabinet.
- White epoxy-coated steel interior insert and exterior casing.
- Single pane energy efficient e-coated glass window door allows inventory observation.
- Easy, hands-free door release.
- Door is fully gasketed and hinged on the right side of the unit.
- Four (4) ¾" non-skid rubber feet.
- Radiant Heat Technology utilizes multiple zone heating that monitors surfaces for safe temperature ranges.
- P-2220 includes a heated center shelf.

Control:




- Electronic control.
- Efficient, balanced heating system allows for low energy consumption and minimal heat loss.
- Adjustable temperature range of 98°F - 180°F (37°C - 82°C).
- Operate in Fahrenheit or Celsius.
- Four digit L.E.D. display.
- On/Off key.
- Actual temperature key.
- Increase and decrease keys.
- Built-in speaker for audible feedback.
- Built-in control button lock-out feature.

Additional features:

- L.E.D. interior lighting casts a comforting blue glow with two (2) different intensity settings and off mode.
- In the event of a power failure, the cabinet will remember its programming and begin to operate as before when power is restored.
- A warming shut-off system, separate from the electronic control, prevents overheating.
- Optional timer key allows the user to program the control to automatically turn on and off at a selected time.



Model	P-2205	P-2210	P-2220
Interior Capacity	1.5 ft³	2.5 ft³	3.5 ft³
Net Weight	56 lbs (25 kg)	67 lbs (30 kg)	91 lbs (41 kg)
Ship Weight	58 lbs (26 kg)	135 lbs (61 kg)	167 lbs (76 kg)
Electrical			
120V Amps	5.0	5.0	6.7
120V kW	0.6	0.6	0.8
Plug (Hospital Grade)	NEMA 5-15P; 15A-125V Plug 	NEMA 5-15P; 15A-125V Plug 	NEMA 5-15P; 15A-125V Plug 
220V Amps	2.7 (Plug - contact factory)	2.7 (Plug - contact factory)	3.6 (Plug - contact factory)
220V kW	0.6 (Plug - contact factory)	0.6 (Plug - contact factory)	0.8 (Plug - contact factory)
230V Amps	2.6	2.6	3.4
230V kW	0.6	0.6	0.8
Plug	BS 1363 Plug (UK only)  	BS 1363 Plug (UK only)  	BS 1363 Plug (UK only)  
	CEE 7/7, 220-230V Plug  	CEE 7/7, 220-230V Plug  	CEE 7/7, 220-230V Plug  

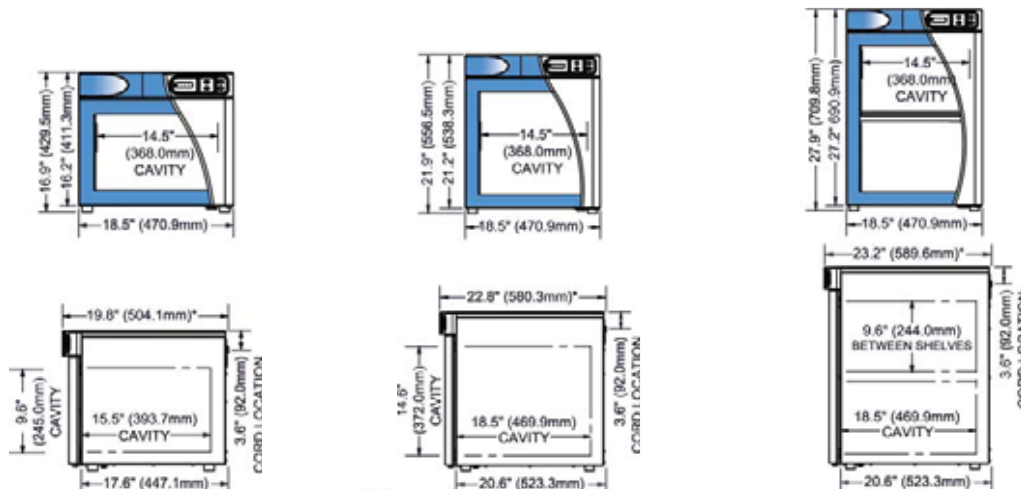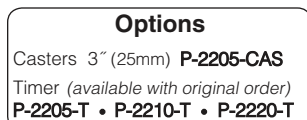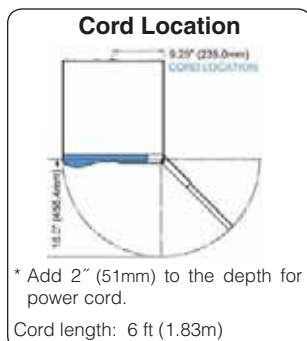



BLANKET & FLUID
WARMING CABINETS

featuring

**RADIANT HEAT
TECHNOLOGY**



FLUID & BLANKET WARMING CABINETS

Deluxe combination blanket & fluid warmers feature 4 digit LED display, window door, control lock programming and auto time programming.

P - 2 1 4 0

FLUID & BLANKET WARMER

COMPARTMENT SIZE

- 7.7 cubic ft. per compartment
- 15.4 cubic ft total
- 18⁷/₈" W x 27" D x 26⁷/₈" H

OVERALL SIZE

- 22¹/₈" W x 34¹⁵/₁₆" D x 72¹/₂" H
- Top compartment for blankets, lower compartment for fluids.
- Fluid compartment has three epoxy-coated baskets each with a maximum capacity of 24 liters. Total capacity is 72 liters.
- Separate temperature controls for each compartment
- Perimeter bumper included

*72 liter
total capacity*

P - 2 1 4 5

FLUID & BLANKET WARMER

COMPARTMENT SIZE

- 38 liter fluid compartment
- 26¹/₂" W x 23³/₄" D x 14⁷/₈" H
- 11.875 cu. ft. blanket compartment
- 26¹/₂" W x 23³/₄" D x 32¹/₂" H

OVERALL SIZE

- 30³/₁₆" W x 27³/₈" D x 70¹⁵/₁₆" H
- Top compartment for fluids, lower compartment for blankets
- Fluid compartment has one easy access epoxy-coated basket with a maximum capacity of 38 1-liter bottles or equivalent
- Separate temperature controls for each compartment

*38 liter
total capacity*



STANDARD MODEL
BLANKET WARMING CABINETS

Standard warmers feature solid door and temperature control.



P - 2 0 1 0 - S
BLANKET WARMER

COMPARTMENT SIZE

- 2.3 cubic feet
- 13 $\frac{3}{4}$ "W x 21 $\frac{1}{2}$ "D x 13 $\frac{7}{8}$ "H

OVERALL SIZE

16 $\frac{3}{16}$ "W x 27 $\frac{1}{2}$ "D x 20 $\frac{5}{8}$ "H

- Analog temperature control
- Ideal for countertop use



P - 2 0 3 0 - S
BLANKET WARMER

COMPARTMENT SIZE

- 7.7 cubic feet
- 17 $\frac{7}{8}$ "W x 26 $\frac{1}{2}$ "D x 26 $\frac{3}{4}$ "H

OVERALL SIZE

22"W x 32 $\frac{1}{2}$ "D x 40 $\frac{5}{8}$ "H

- 3 digit LED display
- Units can be stacked
- 5" casters



P - 2 0 1 0 - T C
TRANSPORT CART

OVERALL SIZE

23 $\frac{3}{16}$ "W x 29"D x 43 $\frac{1}{16}$ "H

- Stainless steel
- Ergonomic Push/Pull handles
- 4" casters with brake and swivel lock
- 300 lbs. total weight capacity at 100 lbs. per shelf
- Compatible with the following warming cabinets:

- **P-2010**
- **P-2010-S**
- **P-2110**



DELUXE MODEL BLANKET WARMING CABINETS

Deluxe blanket warmers feature 4 digit LED display, window door, control lock programming and auto time programming.



P - 2 0 1 0
BLANKET WARMER

COMPARTMENT SIZE

- 2.3 cubic feet
- 13 $\frac{1}{4}$ "W x 21 $\frac{1}{2}$ "D x 13 $\frac{7}{8}$ "H

OVERALL SIZE

16 $\frac{3}{16}$ "W x 25 $\frac{7}{16}$ "D x 20 $\frac{5}{8}$ "H

- Ideal for countertop use
- Unit can be stacked



P - 2 0 2 0
BLANKET WARMER

COMPARTMENT SIZE

- 3.4 cubic feet
- 14 $\frac{1}{4}$ "W x 21 $\frac{5}{8}$ "D x 20 $\frac{1}{8}$ "H

OVERALL SIZE (without bumper)

17 $\frac{5}{8}$ "W x 28"D x 33 $\frac{5}{8}$ "H

- Unit can be stacked
- 5" casters

*Perimeter bumper
optional*



P - 2 0 3 0
BLANKET WARMER

COMPARTMENT SIZE

- 7.7 cubic feet
- 18 $\frac{7}{8}$ "W x 26 $\frac{1}{2}$ "D x 26 $\frac{3}{4}$ "H

OVERALL SIZE (without bumper)

22 $\frac{1}{8}$ "W x 31 $\frac{3}{4}$ "D x 39 $\frac{3}{4}$ "H

- Unit can be stacked
- 5" casters

P - 2 0 4 0
BLANKET WARMER

COMPARTMENT SIZE

- 7.7 cubic feet per chamber (15.4 total cubic feet)
- 18 $\frac{7}{8}$ "W x 26 $\frac{1}{2}$ "D x 26 $\frac{3}{4}$ "H

OVERALL SIZE (without bumper)

22 $\frac{1}{8}$ "W x 31 $\frac{3}{4}$ "D x 73 $\frac{1}{4}$ "H

- Two compartments, each with individual controls
- 5" casters



P - 2 0 5 5
BLANKET WARMER

COMPARTMENT SIZE

- 20.6 cubic feet
- 26 $\frac{1}{2}$ "W x 23 $\frac{3}{4}$ "D x 56 $\frac{3}{4}$ "H

OVERALL SIZE

30 $\frac{1}{4}$ "W x 27 $\frac{1}{4}$ "D x 70 $\frac{7}{8}$ "H

- One compartment with two doors
- 5" casters



FLUID WARMING CABINETS



Deluxe fluid warmers feature 4 digit LED display, window door, control lock programming and auto time programming.

P - 2 1 1 0

FLUID WARMER

COMPARTMENT SIZE

- 2.3 cubic feet
- 13¹/₈"W x 21¹/₂"D x 13¹³/₁₆"H

OVERALL SIZE

16¹/₈"W x 27¹/₂"D x 20¹¹/₁₆"H

- One epoxy-coated basket with a maximum capacity of 16 liters
- Unit can be stacked

16 liter
total capacity



P - 2 1 1 8

FLUID WARMER 18" depth - perfect for built-in OR applications

COMPARTMENT SIZE

- 3.9 cubic feet
- 26¹/₂"W x 17³/₁₆"D x 12¹/₂" H

OVERALL SIZE

30³/₁₆"W x 19³/₈"D x 24⁷/₁₆"H

- One epoxy-coated basket with a maximum capacity of 29 liters
- Two selectable operating modes: 'INJECTION' or 'IRRIGATION'
- Interior illumination



29 liter
total capacity



P - 2 1 2 0

FLUID WARMER

COMPARTMENT SIZE

- 3.4 cubic feet
- 14¹/₈"W x 21⁵/₈"D x 20"H

OVERALL SIZE (perimeter bumper included)

17³/₈"W x 28⁹/₁₆"D x 33⁷/₈"H

- Two epoxy-coated baskets with a maximum capacity of 14 liters per basket
- Unit can be stacked
- 5" casters

28 liter
total capacity



P - 2 1 3 0

FLUID WARMER

COMPARTMENT SIZE

- 7.7 cubic feet
- 18¹⁵/₁₆"W x 26¹/₂"D x 26¹⁵/₁₆"H

OVERALL SIZE (perimeter bumper included)

22"W x 34⁷/₈"D x 40¹¹/₁₆"H

- Three epoxy-coated baskets with a maximum capacity of 24 liters per basket
- Unit can be stacked
- 5" casters

72 liter
total capacity



BLANKET & FLUID WARMERS

COMMON FEATURES

- Electro-thermal Radiant Heat Technology
- Safety coated interior
- Fully insulated
- Optional left hand door swing available
- 125 VAC
- No lint traps
- Heavy duty hardware
- Option of casters or fixed legs

FLUID WARMER FEATURES

- Glass door and perimeter bumper
- 4-digit LED display
- Designed specifically for the warming of irrigation and intravenous fluids
- Individual temperature controls for irrigation fluids and intravenous fluids allow for precise temp settings for the proper fluid type
- Irrigation temperature adjusted from: 98°F to 150°F (37°C to 65.5°C)
Intravenous temperature range: 98°F to 104°F (37°C to 40°C)
- Safety override feature prevents overheating
- Lockout feature



BLANKET WARMER FEATURES

- No fans or other moving parts
- Blankets can be warmed to 200°F (93.3°C)
- Perimeter bumper is optional on Deluxe models
- 4-digit LED display on Deluxe models
- 3 digit LED display on Standard models
- P-2010 models are available with an optional caster stand or leg stand assembly or can be used with P-2010-TC cart
- Safety override feature prevents overheating (*Deluxe models only*)
- Lockout feature (*Deluxe models only*)



RADIANT HEAT TECHNOLOGY

Pedigo warmers feature Radiant Heat Technology to assure accurate set-point temperatures. This patented technology uses an array of electrothermal cables installed in two sides and the bottom of our warmers. Fluid warmers also incorporate the use of a circulation fan to evenly distribute the warmed air to achieve and maintain accurate set-point temperatures. Pedigo fluid warmers have an accuracy of +0/-2°F for injection fluids and an accuracy of +0/-3°F for irrigation fluids. This ensures your fluids are kept at the proper temperature and will always be immediately available at or close to body temperature rather than room temperature without the hook-up time needed for catheter-warming systems.

Fluid Temp Range for **ALL** Fluid Warmers

Intravenous Fluids:	98°F to 104°F	(37°C to 40°C)
Irrigation Fluids:	98°F to 150°F	(37°C to 65.5°C)

Built-in "lock-out" feature prevents inadvertent altering of temperature settings of manufacturer's specifications.

WARRANTY FOR ALL WARMING CABINETS

- One Year Service
- Three Years on Parts
- Lifetime Technical Support



4000 SE Columbia Way • Vancouver, WA 98661



TEL 800.822.3501 • www.pedigo-usa.com



04/10 MM-76



NEW!



WARMWATCH OPTION
FOR FLUID WARMING CABINETS



WARMWATCH OPTION FOR FLUID WARMING CABINETS



Pedigos' injection/irrigation fluid warmers now offer an option with two new features that ensure even greater temperature accuracy and significant labor savings - the WarmWatch option.

Continuous Temperature Verification

Pedigo fluid warmers that are equipped with the optional WarmWatch package feature a continuous temperature monitoring system that provides users instant visual confirmation of the chamber temperature. In addition to the standard control probe, an independent monitoring probe located very close to the control probe verifies that the chamber reaches and maintains the set-point temperature with $\pm 0.1^\circ\text{C}$ ($\pm 0.2^\circ\text{F}$) accuracy from $37 - 40^\circ\text{C}$ ($98 - 104^\circ\text{F}$) and $\pm 0.17^\circ\text{C}$ ($\pm 0.3^\circ\text{F}$) accuracy above $40 - 66^\circ\text{C}$ ($104 - 150^\circ\text{F}$). Two independent LED displays on the control panel show the set-point and actual interior temperature at a glance.


(The chambers are designed not to exceed the set-point temperature, thus the nominal actual temperature will always be one to one-and-a-half degrees lower than the set-point temperature.)

Digital Documentation

In addition to continuous temperature monitoring, the WarmWatch option provides digital documentation of the warmer's historical temperature performance data. This offers significant advantages over manually recording temperature set-points or control probe displays since these values are not independently obtained, and do not verify the actual chamber temperature. The independent monitoring probe in this option ensures that the actual chamber temperature is recorded. Automatic collection also ensures consistent data samples.

This recorded data can be downloaded to a USB flash drive and then to a computer at any desired time interval—up to six months between downloads. Normal warming functions continue to operate while the data is downloading. The data download to the USB

drive includes both historical temperature information and a detailed performance summary for the warmer. These reports provide convenient and reliable data analysis for the life of the warmer.

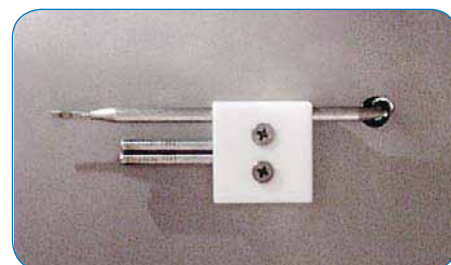

Hi-limit sensor (top), independent monitoring probe (middle), and control probe (bottom)



Digital Documentation Features

- Stores up to six months of data and displays a warning when the unit approaches its memory limit
- Downloading begins automatically when a USB flash drive is connected to unit
- Accuracy report includes:
 - ✓ Warmer model number, serial number, and date of manufacture
 - ✓ Date of report/download
 - ✓ Software version
 - ✓ Frequency of readings used in evaluation
 - ✓ Dates included in evaluation
 - ✓ Number of days included in report
 - ✓ Number of readings evaluated for accuracy calculation
 - ✓ Number of readings excluded from accuracy calculation due to a door opening within two hours of reading
 - ✓ Mode of fluid warming (injection or irrigation)
 - ✓ Set-point temperature
 - ✓ Average temperature
 - ✓ Temperature range during evaluation time period
 - ✓ Accuracy evaluation indicating whether unit successfully maintained accuracy parameters while set to selected set point
- Raw temperature data reported in a comma-delimited (CSV) file

Continuous Temperature Verification Features

- Temperature monitoring via an independent monitoring probe
- Displays set-point and actual temperature simultaneously



STRETCHER



P - 1 7 0

OVERALL SIZE

78 $\frac{1}{4}$ " L x 28" W x 33" H

FEATURES

- Designed as a cost efficient transport stretcher
- Heavy-gauge, 1 $\frac{3}{8}$ " O.D. tubular frame
- Continuous heavy rubber bumper
- Convenient hand holds on each corner
- Long-life polyurethane, double ball bearing, 8" swivel casters for quiet movement



P - 1 7 1

OVERALL SIZE

78 $\frac{1}{4}$ " L x 28" W x 33" H

FEATURES

- Designed as a cost efficient transport stretcher
- Two combination swivel lock and brake lock casters mounted on diagonal corners for easier control
- Heavy-gauge, 1 $\frac{3}{8}$ " O.D. tubular frame
- Continuous heavy rubber bumper
- Convenient hand holds on each corner
- Long-life polyurethane, double ball bearing, 8" swivel casters for quiet movement

OPTIONAL ACCESSORIES

- 2 $\frac{1}{2}$ " thick black mattress pad with soft vinyl cover
- Combination swivel lock & brake lock casters
- O₂ holder for E-size tank (*requires lower shelf*)
- IV pole (*includes holder and storage tube*)
- Stainless steel fowler head rest (0° to 45° range)
- 7 $\frac{1}{2}$ " rolling rubber bumpers
- Restraint straps
- Lower shelf
- Side rails



Stainless Steel Fowler Head Rest




7 $\frac{1}{2}$ " Rolling Bumpers



Swivel Lock & Brake Lock Caster




Restraint Straps

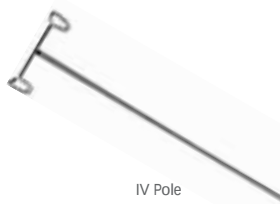

IV Pole



2 $\frac{1}{2}$ " Thick Mattress Pad

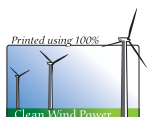



500

Pediatric Crib

Stretcher



PEDIATRIC STRETCHER/CRIB

Function as Appealing as its Form

For the care of your young patients, Pedigo's 500 Pediatric Crib Stretcher with its full function capabilities represents the standard of excellence in hospital cribs. The aesthetically pleasing design makes it perfect for any facility.

As a pediatric transport stretcher, the 500 rolls easily on 5" oversized casters. Two of the four casters lock to secure the stretcher while a third steering caster may be engaged for increased control and maneuverability.

Optional IV poles allow for bottles or bags to be suspended inside or outside the stretcher area.

The unique rounded corners reduce the overall size of the 500, yet provide the same amount of useful area on the 26½" x 55¾" mattress. An optional

utility shelf located at the base may be used for storing blankets or the child's personal items.

Total Accessibility

The two-post support system provides 97% access to the patient by the medical team.

Fixed Security Top

With the optional translucent security top in position, the 500 becomes a full function pediatric crib. The molded top has an exclusive quick disconnect for fast, easy removal.



Optional Departmental Color Coding

Identify your 500s with the departmental color of your choice on each end.

A True One-Handed Crib

Dual-action release handle is childproof and easy to operate, even with a child in one arm. Side rails have two-inch bar spacing and a special epoxy finish that is warmer to the touch than chrome.



For complete containment of the child, simply raise the counter balanced side rails to enclose the crib. Because the rails are counter balanced, there will be no sudden drop when released, preventing injury to the child or provider.

500 STANDARD FEATURES

NEW!

- Two 5" locking casters and one steering caster for increased maneuverability
- Trendelenburg/Reverse Trendelenburg
- One handed rail operation
- White powder coat finish
- Unique rounded corners
- Rails with 2" bar spacing
- Two IV pole receptacles
- Full function fowler
- Litter Top: Solid pan
- Color decals: Blue

✓ LATEX-FREE



SPECIFICATIONS

Model Number	500 Stretcher Crib
Patient Surface	26½" x 55¾" (67.3 cm x 141.6 cm)
Overall Length	60" (152.4 cm)
Overall Width	30¼" (76.8 cm)
Height	
Fixed Litter Top	32¼" (81.9 cm)
Top of Crib	59½" (151.1 cm)
With Security Top	74" (188 cm)
Litter Positioning	
Fowler Angles	0°, 20°, 30°
Trendelenburg	0°, 5°, 10°
Reverse Trendelenburg	-5°
Side Rail Intermediate Stop	13¾ (39.9 cm)
Casters	5" (12.7 cm)
Warranty	One year limited

OPTIONS

- | | |
|--|--|
| 5103001 Utility shelf | 5114001 Oxygen tank holder |
| 5104002 Mounting kit for traction adapter and croupette adapter | 5116001 IV Pole |
| 5106001 Infusion pump adapter | 5116001 Mattress: Foam pad with Staph CHEK™ cover |
| 5107001 Windowed Security top with 6 mounting clips (will accommodate a mist tent) | 5117001 Security top, translucent, with 6 mounting clips (will accommodate a mist tent) |
| 5108001 Croupette humidifier adapter (requires 5104002 for mounting) | 5119001 Traction adapter, pair (includes mounting kit) |
| 5109002 Crib bumper pad | 5817002 Pressure care pad, 4" |
| | Color decals Red Green Blue |





stretcher series

5110

general transport stretcher



stretcher series

5110

The full function 5110 Series general transport stretcher provides a versatile yet economical patient transport solution. This is a non-hydraulic adjustable stretcher.



1

RAIL SYSTEM:

Easy to wipe clean and disinfect, full length tuck-away side rails allow for easier patient transfer by creating a zero transfer gap. Rails feature a release latch centrally located on either side of the litter top under the patient surface.



2

PNEUMATIC ASSIST FOWLER:

Provides additional assistance for the safe and easy raising and lowering of the patient.



3

TRENDELENBURG:

A single pneumatic cylinder assists in placing the litter top in Trendelenburg/ reverse Trendelenburg position.



4

CENTER COLUMN SUPPORT

Sturdy steel column design with five fixed litter top height positions to choose from: 28", 29", 30", 31" or 32".



5

TRUE DIRECTION STEERING (TDS)

For improved stretcher control engage the True Direction Steering system. The steer caster located at the patient's right foot end of the stretcher will lock once the steering pedal is applied.

stretcher options




FOLD DOWN IV POLE



Available as an option, the fold down IV pole tucks out of the way when not needed and remains attached to the stretcher when in use.



IV CADDY

Quickly and easily secure an IV pole to a stretcher using Pedigo's IV Caddy. Allows one caregiver to manage the safe transport of the patient, stretcher and IV equipment. Stores under litter top when not in use.

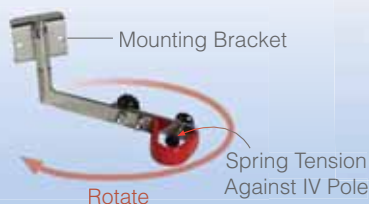

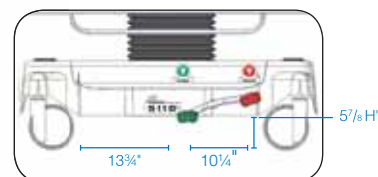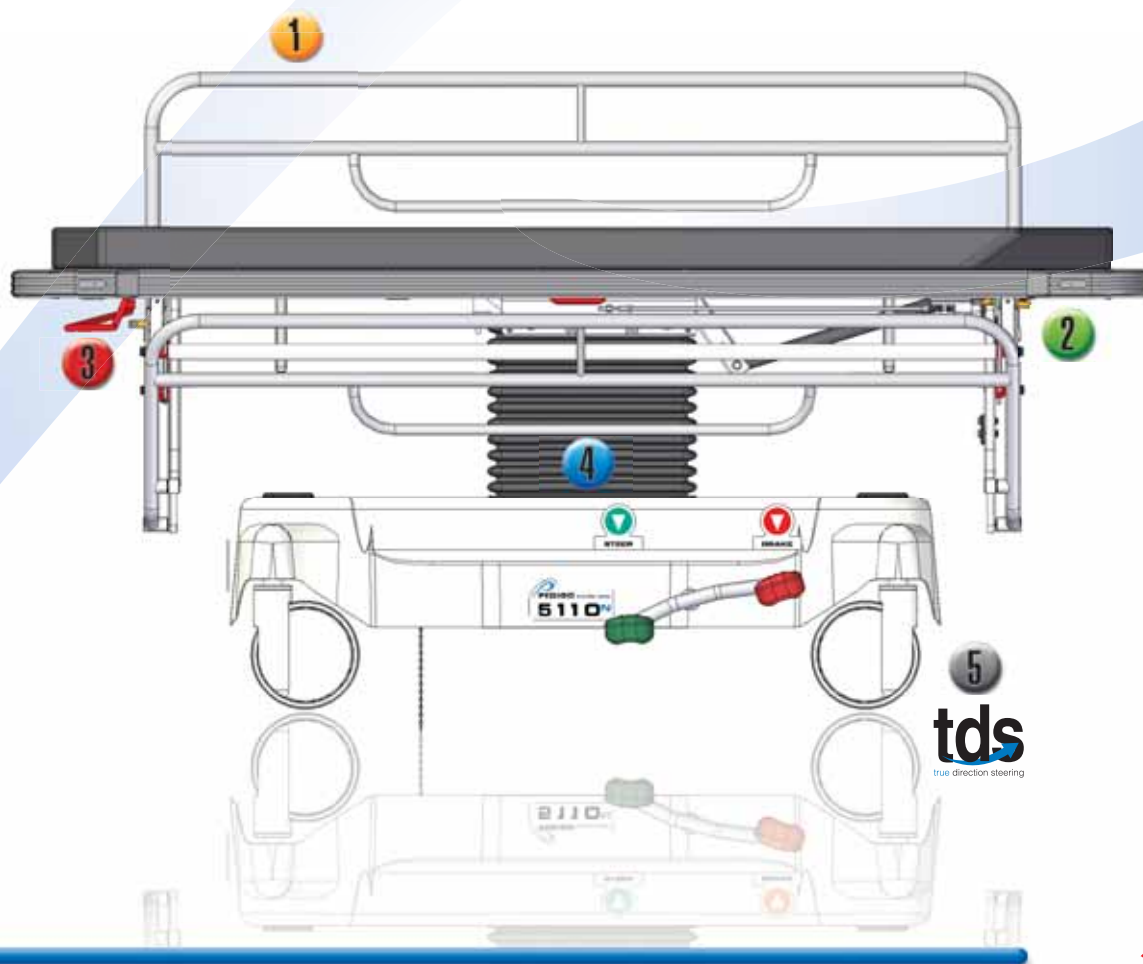

PATIENT LIFTS:

The 5110 Series Stretcher has 57 $\frac{7}{8}$ " H" of clearance between floor and cowl to easily accommodate and utilize patient lifts.

✓ **LATEX-FREE**



ZERO TRANSFER GAP

Tuck-Away side rails allows the caregiver to position the 5110 Series Stretchers flush with the bed for a zero transfer gap and easier patient transfer. No expensive or cumbersome transfer devices are required.



CuVerro®
Antimicrobial Copper Surfaces

Now as an added protection in the ongoing battle with bacteria, Pedigo offers the option of side rails made from **CuVerro®** antimicrobial copper alloy.

Copper is a natural anti-microbial element and **CuVerro®** is the ONLY solid touch surface material registered by the EPA to continuously kill bacteria¹ that cause infections and pose a risk to human health. Lab tests show 99.9% of bacteria¹ associated with disease and infection, including MRSA¹, are killed within two hours when in contact with copper. Bacteria will continue to be killed 24 hours a day, week after week, between regular cleanings.

¹ CuVerro® alloys have proven effectiveness on MRSA (Methicillin-Resistant Staphylococcus aureus), Staphylococcus aureus, Enterobacter aerogenes, Pseudomonas aeruginosa, and Escherichia coli 0157:H7. CuVerro is a registered trademark of GBC Metals LLC, d/b/a Olin Brass.

SPECIFICATIONS

Model Number	5110 ^N	5110 ^W
Weight Capacity	500 lbs. (227 kg)	500 lbs. (227 kg)
Overall Length	83" (211 cm)	83" (211 cm)
Overall Width	29.5" (75 cm)	33.5" (85 cm)
Fixed Height Range		
Low	28" (71 cm)	28" (71 cm)
	29" (73 cm)	29" (73 cm)
	30" (76 cm)	30" (76 cm)
	31" (79 cm)	31" (79 cm)
High	32" (81 cm)	32" (81 cm)
Litter Positioning		
Fowler Backrest	0°-90°	0°-90°
Trendelenburg	+12°	+12°
Reverse Trendelenburg	-7°	-7°
Patient Surface	77" x 27" (196 x 69 cm)	77" x 31" (196 x 77 cm)
Side Rails	15" x 63 ⁵ / ₈ " (38 x 162 cm)	15" x 63 ⁵ / ₈ " (38 x 162 cm)
Casters	8" (20 cm)	8" (20 cm)
Warranty	1 year parts and labor	1 year parts and labor

STANDARD FEATURES

- 500 lbs. weight capacity
- True Direction Steering (TDS)
- Side rails - full length, single release
- Center column patient support
- Non-hydraulic
- 5 fixed litter top height options: 28" to 32"
- A single pneumatic cylinder assists in placing the litter top in Trendelenburg/ reverse Trendelenburg position.
- Pneumatic assist fowler 400 N (90#) force
- Brake & Steer activated foot pedals
- Four corner IV pole receptacles
- Four restraint strap positions
- Full PVC perimeter bumper
- 8" easy mobility casters
- Latex-free

OPTIONS

Accessory	N	W	Accessory	N	W	Accessory	N	W
X-Ray Board			IV Pole Options:			Additional Options:		
Full	750203	750204	Fixed Height	2101	same	Bed Exit Monitor	750216	same
Partial	750205	750206	Adjustable / Removable ..	5984001	same	Storage Bag Hook	750215	same
Mattresses:			Adjustable / Attached	750501	750502	Paper Roll Holder	5966001	5966003
3" Standard (black)	5828002	5831002	Fold Down	59175001	same	Quick Release O ₂ Holder ..	750213	same
3" Premium (blue)	5850002	5852002	IV Pole Caddy	59163001	same	O ₂ Tank Holder (Vertical) ..	5981001	same
4" Standard (black)	5829002	5832002	Push/Pull Bar Options:			Head/Foot Extension	5996001	same
4" Premium (blue)	5851002	5853002	Push/Pull Bar (Removable) ..	5985001	5986001	Head/Foot 3" Pad	5708001	same
4" Pressure Relief	5833002	5834002	Bar Pad	5704001	5705001	Stirrups Set	59100001	same
5" Premium (blue)†	5855002	5856002	Chart Holder	5967001	same	Color-Coded Bumper Options:		
† Special order only.			Monitor Tray Board	59115001	59116001	Red	750401	same
Restraint Straps Options:			Knee Flex			Teal	750402	same
Velcro (pair)	5903001	same	Hand Crank	750201	750202			
Buckle (pair)	5904001	same	Hydraulic	750208	750209			
Side Rail Options:								
CuVerro® Rails	750303	same				Red		
Rail Pad (pair)	5716001	same				Teal		
			Leg Lift			Grey (standard)		
			Knee Flex					
			Trauma					





stretcher series
5400
universal procedure stretcher



stretcher series 5400

The 5400 Series universal transport stretcher is versatile, yet economical. It's the right stretcher for ED, OR, PACU, or any department where cost, efficiency and quality are required.



1

RAIL SYSTEM:

Easy to wipe clean and disinfect, full length tuck-away side rails allow for easier patient transfer by creating a zero transfer gap. Rails feature a release latch centrally located on either side of the litter top under the patient surface.



2

PNEUMATIC ASSIST FOWLER:

Provides additional assistance for the safe and easy raising and lowering of the patient.



3

HYDRAULIC HEIGHT ADJUSTMENT FOOT PEDALS:

Hydraulic Foot Pedals are located on each side of the stretcher within the patient care zone. Color-coded pedals and function graphics are quickly and easily understood.



4

CENTER COLUMN SUPPORT:

Sturdy steel column with single hydraulic lift cylinder design uses sixteen precision ball-bearings to smoothly raise and lower the patient.



5

TRUE DIRECTION STEERING (TDS)

For improved stretcher control engage the True Direction Steering system. The steer caster located at the patient's right foot end of the stretcher will lock once the steering pedal is applied.

stretcher options




FOLD DOWN IV POLE



Available as an option, the fold down IV pole tucks out of the way when not needed and remains attached to the stretcher when in use.



IV CADDY

Quickly and easily secure an IV pole to a stretcher using Pedigo's IV Caddy. Allows one caregiver to manage the safe transport of the patient, stretcher and IV equipment. Stores under litter top when not in use.

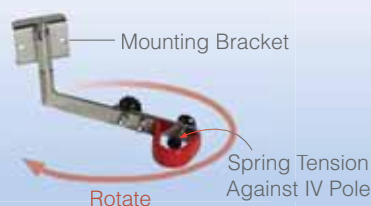

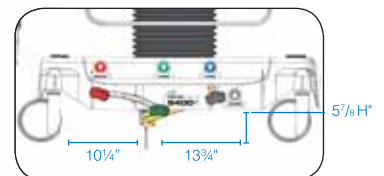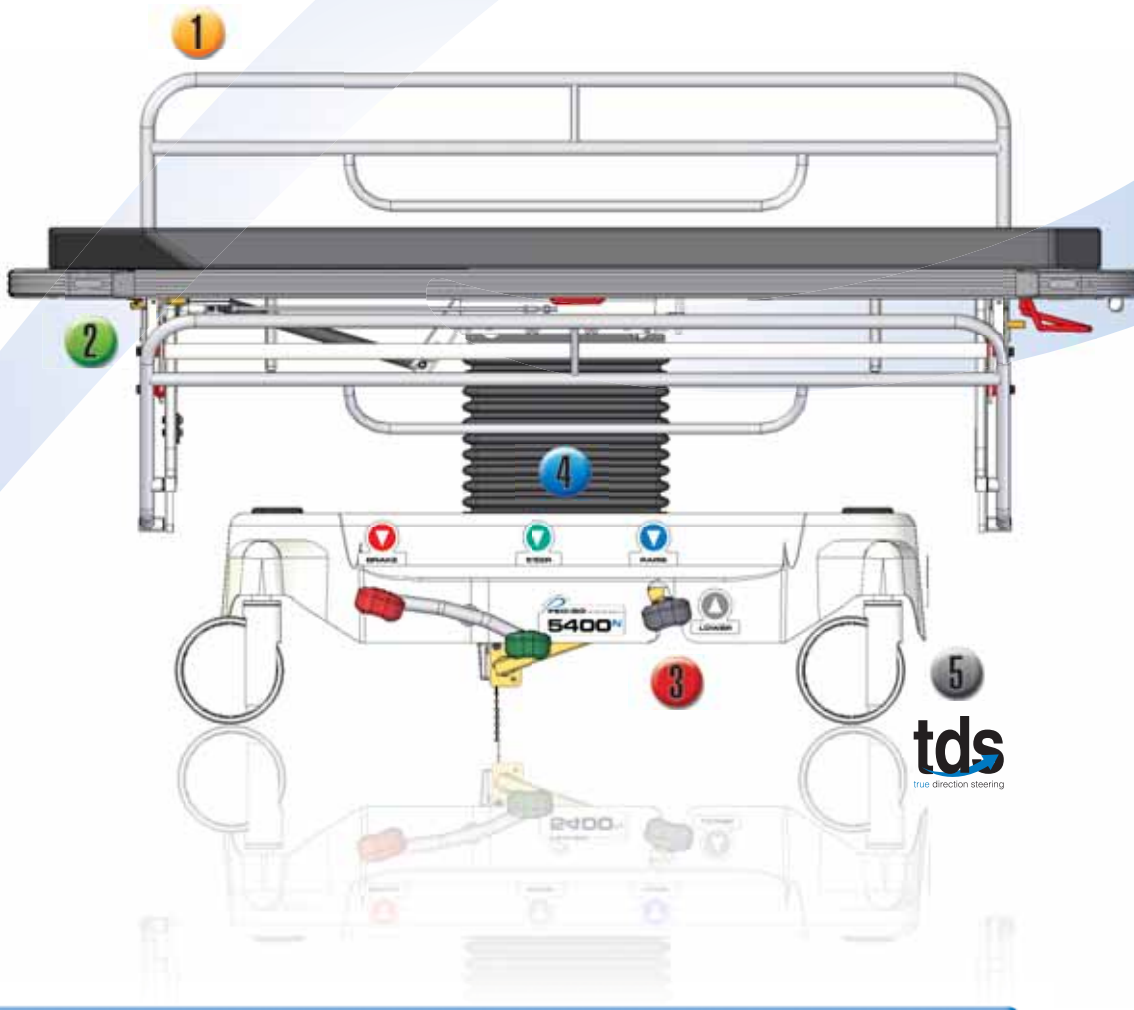

PATIENT LIFTS:

The 5400 Series Stretcher has 5⁷/₈ H" of clearance between floor and cowling to easily accommodate and utilize patient lifts.

✓ **LATEX-FREE**



ZERO TRANSFER GAP

Tuck-Away side rails allows the caregiver to position the 5400 Series Stretchers flush with the bed for a zero transfer gap and easier patient transfer. No expensive or cumbersome transfer devices are required.

CuVerro®
Antimicrobial Copper Surfaces

Now as an added protection in the ongoing battle with bacteria, Pedigo offers the option of side rails made from **CuVerro®** antimicrobial copper alloy.

Copper is a natural anti-microbial element and **CuVerro®** is the ONLY solid touch surface material registered by the EPA to continuously kill bacteria¹ that cause infections and pose a risk to human health. Lab tests show 99.9% of bacteria¹ associated with disease and infection, including MRSA¹, are killed within two hours when in contact with copper. Bacteria will continue to be killed 24 hours a day, week after week, between regular cleanings.

¹ CuVerro® alloys have proven effectiveness on MRSA (Methicillin-Resistant Staphylococcus aureus), Staphylococcus aureus, Enterobacter aerogenes, Pseudomonas aeruginosa, and Escherichia coli 0157:H7. CuVerro is a registered trademark of GBC Metals LLC, d/b/a Olin Brass.



SPECIFICATIONS

Model Number	5400 ^N	5400 ^W
Weight Capacity	500 lbs. (227 kg)	500 lbs. (227 kg)
Overall Length	83" (211 cm)	83" (211 cm)
Overall Width	29.5" (75 cm)	33.5" (85 cm)
Height Range		
Low	23" (58 cm)	23" (58 cm)
High	32.75" (83 cm)	32.75" (83 cm)
Litter Positioning		
Fowler Backrest	0°-90°	0°-90°
Trendelenburg	+12°	+12°
Reverse Trendelenburg	-7°	-7°
Patient Surface	77" x 27" (196 x 69 cm)	77" x 31" (196 x 77 cm)
Side Rails	15" x 63 ⁵ / ₈ " (38 x 162 cm)	15" x 63 ⁵ / ₈ " (38 x 162 cm)
Casters	8" (20 cm)	8" (20 cm)
Warranty	1 year parts and labor	1 year parts and labor

STANDARD FEATURES

- 500 lbs. weight capacity
- True Direction Steering (TDS)
- Side rails - full length, single release
- Center column design with precision ball bearings
- Positive/negative Trendelenburg single pneumatic assist cylinder
- Pneumatic assist fowler 400 N (90#) force
- Dual side-mounted hydraulic foot pedals
- Brake & Steer activated foot pedals
- Four corner IV pole receptacles
- Four restraint strap positions
- Full PVC perimeter bumper
- 8" easy mobility casters
- Latex-free

OPTIONS

Accessory	N	W	Accessory	N	W	Accessory	N	W
X-Ray Board			IV Pole Options:			Additional Options:		
Full	750203	750204	Fixed Height	2101	same	Digital Scale	750214	same
Partial	750205	750206	Adjustable / Removable	5984001	same	Bed Exit Monitor	750216	same
Mattresses:			Adjustable / Attached	750501	750502	Extra Release Pedal	750105	same
3" Standard (black)	5828002	5831002	Fold Down	59175001	same	(TDS True Direction Steering)		
3" Premium (blue)	5850002	5852002	IV Pole Caddy	59163001	same	iS6™ Instant Steer 6 Wheel	750103	same
4" Standard (black)	5829002	5832002	Push/Pull Bar Options:			Extra Release Pedal	750106	same
4" Premium (blue)	5851002	5853002	Push/Pull Bar (Removable)	5985001	5986001	(iS6 Six Wheel Steering)		
4" Pressure Relief	5833002	5834002	Bar Pad	5704001	5705001	Storage Bag Hook	750215	same
5" Premium (blue)†	5855002	5856002	Chart Holder	5967001	same	Paper Roll Holder	5966001	5966003
† Special order only.			Monitor Tray Board 59115001 59116001			Quick Release O ₂ Holder	750213	same
Restraint Straps Options:			(includes chart holder 5967001)			O ₂ Tank Holder (Vertical)	5981001	same
Velcro (pair)	5903001	same	Knee Flex			Head/Foot Extension	5996001	same
Buckle (pair)	5904001	same	Hand Crank	750201	750202	Head/Foot 3" Pad	5708001	same
Side Rail Options:			Hydraulic	750208	750209	Stirrups Set	59100001	same
CuVerro® Rails	750303	same				Color-Coded Bumper Options:		
Rail Pad (pair)	5716001	same				Red	750401	same
						Teal	750402	same
								
								
								
						Red	Teal	Grey (standard)





stretcher series
7500
transport / trauma stretcher



iS6™
six wheel steering

stretcher series 7500

Pedigo's 7500 Series Stretcher allows for the transfer of patients easily and comfortably with a load capacity up to 750 pounds. Superior stability, mobility, and braking capabilities make the 7500 versatile for all your department needs such as OR, ED, PACU, Patient Transportation, Radiology or Cath Lab.



1

RAIL SYSTEM:

Full length stainless steel side rails feature integrated push-pull handles with latch releases located at all four corners of the stretcher. Easy one-handed operation. No hinges, mechanical connections or plastic couplings to trap debris or fluid. Easy to wipe clean and disinfect.



2

OXYGEN TANK HOLDER:

Pedigo's ergonomic Quick-Release O₂ Tank Holder is side mounted and conveniently located under the patient surface to reduce caregiver bending. (US Patent 7,865,984 B2)



3

new

IS6™ - INSTANT STEER SIX WHEEL STEERING:

The two IS6™ wheels are centered directly under the chassis for increased maneuverability and control. Dual ball bearing wheels feature a softer wheel tread for improved grip and traction.



4

CENTER COLUMN SUPPORT:

Sturdy steel column with single hydraulic lift cylinder design uses sixteen precision ball-bearings to smoothly raise and lower the patient.



5

HYDRAULIC HEIGHT ADJUSTMENT FOOT PEDALS:

Hydraulic Foot Pedals are located on each side of the stretcher within the patient care zone. Color-coded pedals and function graphics are quickly and easily understood.



STRETCHER SCALE

Back lit LCD scale display reads to the nearest pound with a +/- margin of 5 pounds throughout the entire 750 pound range. Instrumentation cable is routed through frame so as not to interfere with patient or peripheral components.



ZERO

Scale can be zeroed to adjust for any accessories or additional gear on the litter top.



WEIGH

Stores patient weight to memory when patient is weighed.



RECALL /
HOLD

Weight can be recalled even if patient is not on stretcher.





MODE

One button switch from pounds to kilograms.

stretcher options



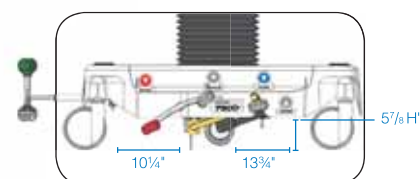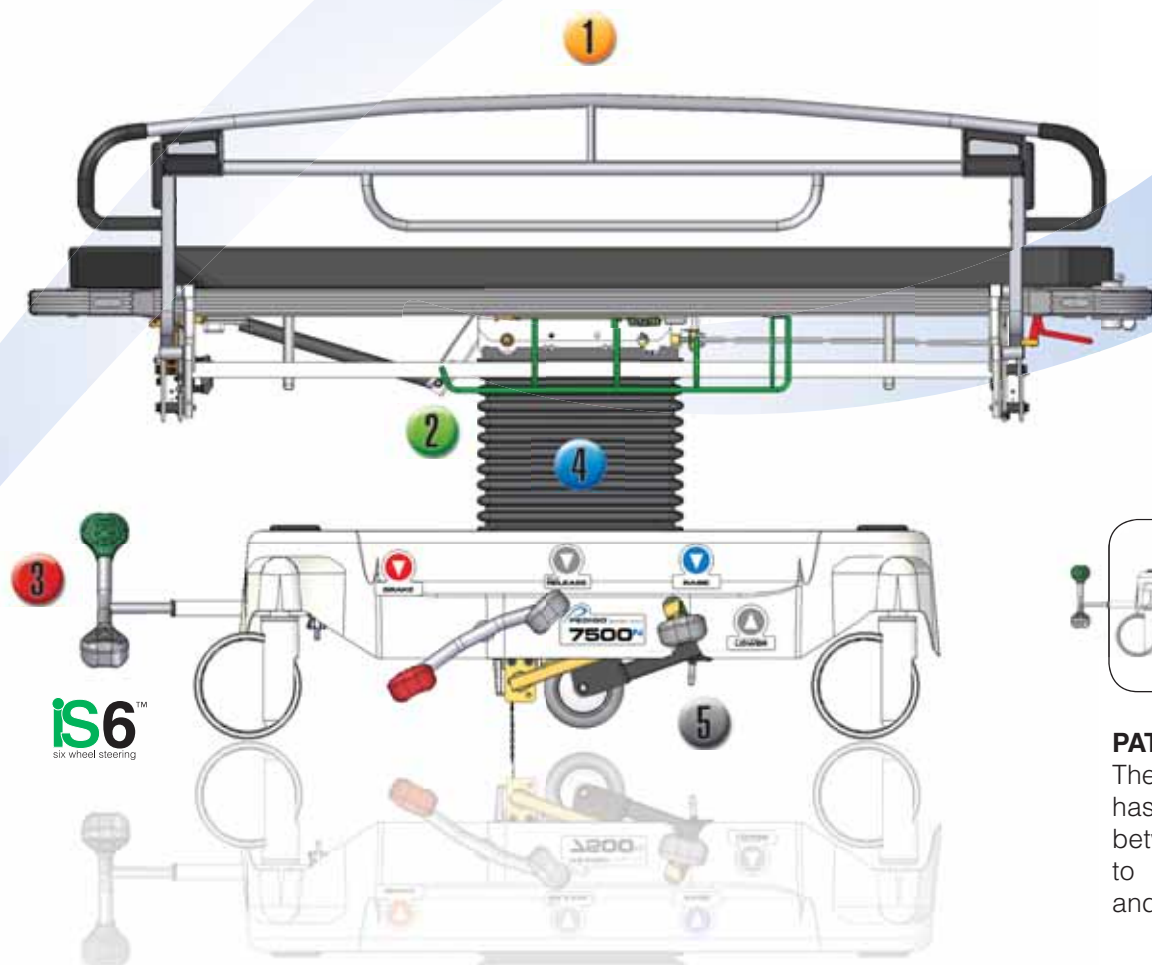

PATIENT LIFTS:

The 7500 Series Stretcher has 5⁷/₈ H" of clearance between floor and cowl to easily accommodate and utilize patient lifts.

✓ LATEX-FREE

Antimicrobial Copper **Cu⁺**

CuVerro®

Antimicrobial Copper Surfaces

Now as an added protection in the ongoing battle with bacteria, Pedigo offers the option of side rails made from **CuVerro®** antimicrobial copper alloy.

Copper is a natural anti-microbial element and **CuVerro®** is the ONLY solid touch surface material registered by the EPA to continuously kill bacteria¹ that cause infections and pose a risk to human health. Lab tests show 99.9% of bacteria¹ associated with disease and infection, including MRSA¹, are killed within two hours when in contact with copper. Bacteria will continue to be killed 24 hours a day, week after week, between regular cleanings.

¹ CuVerro® alloys have proven effectiveness on MRSA (Methicillin-Resistant Staphylococcus aureus), Staphylococcus aureus, Enterobacter aerogenes, Pseudomonas aeruginosa, and Escherichia coli 0157:H7. CuVerro is a registered trademark of GBC Metals LLC, d/b/a Olin Brass.



IV CADDY

Quickly and easily secure an IV pole to a stretcher using Pedigo's IV Caddy. Allows one caregiver to manage the safe transport of the patient, stretcher and IV equipment. Stores under litter top when not in use.



SPECIFICATIONS

Model Number	7500 ^N	7500 ^W
Weight Capacity	750 lbs. (340 kg)	750 lbs. (340 kg)
Overall Length	83" (211 cm)	83" (211 cm)
Overall Width	29.5" (75 cm)	33.5" (85 cm)
Height Range		
Low	23" (58 cm)	23" (58 cm)
High	32.75" (83 cm)	32.75" (83 cm)
Litter Positioning		
Fowler Backrest	0°-90°	0°-90°
Trendelenburg	+12°	+12°
Reverse Trendelenburg	-7°	-7°
Patient Surface	77" x 27" (196 x 69 cm)	77" x 31" (196 x 77 cm)
Side Rails	14" x 76" (36 x 193 cm)	14" x 76" (36 x 193 cm)
Casters	8" (20 cm)	8" (20 cm)
Warranty *	1 year parts and labor	1 year parts and labor

* 7500 Special Package Warranty (7500-SPEC): 2 year parts and labor.

STANDARD FEATURES

- 750 lbs. weight capacity
- **IS6™** Instant Steer 6 Wheel Steering
- Side Rails - dual release stainless steel glide & tuck
- Integrated push/pull handles
- Center column design with precision ball bearings
- Quick-Release O₂ Tank Holder
(US Patent 7,865,984 B2)
- Positive/negative Trendelenburg dual pneumatic assist cylinders
- Mechanical Level Lock
- Pneumatic assist fowler 600N (135#) force
- Dual side-mounted hydraulic foot pedals
- Brake & Steer activated foot pedals
- Four corner IV pole receptacles
- Four restraint strap positions
- Full PVC perimeter bumper
- 8" easy mobility casters
- Latex-free

OPTIONS

Accessory	N	W	Accessory	N	W	Accessory	N	W
X-Ray Board			IV Pole Options:			Additional Options:		
Full	750203	750204	Fixed Height	2101	same	Digital Scale	750214	same
Partial	750205	750206	Adjustable / Removable	5984001	same	Extra Release Pedal	750106	same
Mattresses:			Adjustable / Attached	750501	750502	(IS6 Six Wheel Steering)		
3" Standard (black)	5828002	5831002	Fold Down	59175001	same	Storage Bag Hook	750215	same
3" Premium (blue)	5850002	5852002	IV Pole Caddy	59163001	same	Paper Roll Holder	5966001	5966003
4" Standard (black)	5829002	5832002	Push/Pull Bar Options:			O ₂ Tank Holder (Vertical) ..	5981001	same
4" Premium (blue)	5851002	5853002	Push/Pull Bar (Removable) ..	5985001	5986001	400N (90#) pneumatic		
4" Pressure Relief	5833002	5834002	Bar Pad	5704001	5705001	assist fowler	750602	same
5" Premium (blue)†	5855002	5856002	Chart Holder	5967001	same	Head/Foot Extension	5996001	same
† Special order only.			Monitor Tray Board	59115001	59116001	Head/Foot 3" Pad	5708001	same
Restraint Straps Options:			(includes chart holder 5967001)			Stirrups Set	59100001	same
Velcro (pair)	5903001	same	Knee Flex			Color-Coded Bumper Options:		
Buckle (pair)	5904001	same	Hand Crank	750201	750202	Red	750401	same
Side Rail Options:			Hydraulic	750208	750209	Teal	750402	same
CuVerro® Rails	750302	same	  			  		
Rail Pad (pair)	5716001	same	Leg Lift Knee Flex Trauma			Red Teal Grey (standard)		





- Copper is a natural anti-microbial element and the only solid touch surface material registered by the EPA to continuously kill bacteria¹ that cause infections and pose a risk to human health.
- In the CuVerro® alloys registered with the EPA, copper is the active bacteria¹-killing ingredient. Lab tests show 99.9% of bacteria¹ associated with disease and infection, including MRSA¹, are killed within two hours when in contact with copper. Bacteria will continue to be killed 24 hours a day, week after week, between regular cleanings.
- Beyond EPA claims, clinical trials and research have shown copper alloys to be effective in a wider range of properties:
 Additional tests have shown a “halo” effect with the alloy. i.e. non- copper surfaces such as wood adjacent to CuVerro® surfaces, have shown up to an 87% reduction in bio-burden.²
 Studies on copper alloys containing greater than 70% copper demonstrate significant reduction in survival of C.diff. (Clostridium difficile) spores versus that of stainless steel. Complete death of spores occurred within 24 to 48 hours versus no significant death rate on stainless steel after 168 hours.³
- Other than Copper alloys, no other solid touch surface materials have been registered by the EPA for their antimicrobial properties. Copper is the only material that is EPA-registered as an anti-microbial solid surface material.
- Unlike coatings, the antimicrobial effectiveness of CuVerro® will not wear off over time. The antimicrobial properties of CuVerro® surfaces will not degrade from scuffs, scratches, regular cleaning or recommended use.
- Leading hospitals are involved in bio-burden studies funded by the DOD that are showing reduction in bioburden levels that approach reductions seen in laboratory settings. Pedigo will be participating future clinical trials.
- Studies monitoring bacteria¹ load on ICU surfaces in regularly cleaned rooms have shown that bed rails have the highest levels of bacteria¹ of all the surfaces sampled.⁴ Rails made with CuVerro® significantly reduce risks associated with this bioload.
- Pedigo rails are designed to be easily and thoroughly cleaned. CuVerro® antimicrobial surfaces continue to kill bacteria¹ between normal cleanings offering additional protection.

NOTE:

The use of copper alloy surfaces is a supplement to and not a substitute for standard infection control practices; users must continue to follow all current infection control practices, including those practices related to cleaning and disinfection of environmental surfaces.

¹ CuVerro® alloys have proven effectiveness on MRSA (Methicillin-Resistant Staphylococcus aureus), Staphylococcus aureus, Enterobacter aerogenes, Pseudomonas aeruginosa, and Escherichia coli O157:H7. CuVerro is a registered trademark of GBC Metals LLC, d/b/a Olin Brass.

² Copper Surfaces Reduce Microbial Burden in Out-Patient Infectious Disease Practice

B E Hirsch, H Attaway, R Nadan, S Fairey, J Hardy, G Miller, S Rai, D Armellino, M Schilling, W Moran, P Sharpe, A Estelle, J H Michel, H T Michels, M G Schmidt - presented in a poster session at the Interscience Conference on Antimicrobial Agents and Chemotherapy (ICAAC) in Boston, MA, September 13 2010.

³ Antimicrobial Efficacy of Copper Surfaces against Spores and Vegetative Cells of Clostridium difficile: the germination theory L. J. Wheeldon, T. Worthington, P. A. Lambert, A. C. Hilton, C. J. Lowden and T. S. J. Elliott, Journal of Antimicrobial Chemotherapy 2008 62(3):522-525; doi:10.1093/jac/dkn219.

⁴ Microbial Burden of Objects in ICU Rooms C D Salgado, K A Sepkowitz, T Plaskett, J F John, J R Cantey, H H Attaway, L L Steed, H T Michels, M G Schmidt. October 2008.

⁵ Effectiveness of Copper Contact Surfaces in Reducing the Microbial Burden (MB) in the Intensive Care Unit (ICU) of Hospital del Cobre, Calama, Chile V Prado, C Durán, M Cresto, A Gutierrez, P Sapiain, G Flores, H Fabres, C Tardito, M Schmidt. Poster 56.044, presented at the 14th International Conference on Infectious Diseases, Miami, March 11, 2010.



750NE

Series Stretcher



STANDARD FEATURES

- 750 lbs. weight capacity rating
- Quick-Release O₂ Tank Holder (US Patent 7,865,984 B2)
- Anodized aluminum glide side rails
- Integrated push/pull handles with latch release
- Dual side-mounted hydraulic raise/lower foot pedals
- Single column hydraulic cylinder design
- Pneumatic assist positive and negative Trendelenburg
- Pneumatic assist Fowler backrest
- Knee Flex/Leg Lift - hydraulic
- Brake & Steer foot pedals
- Four corner IV pole receptacles
- Four restraint strap positions
- Instant Steer
- Head Rest Positioning
- 4" Premium mattress
- 8" easy mobility casters
- Latex-free

OPTIONS

Factory Installed Options:

IV Pole Options:

- Fixed
- Adjustable / Removable
- Adjustable / Attached
- Fold Down
- IV Pole Caddy

Restraint Straps Options:

- Velcro
- Buckle

Fowler Options:

- Pneumatic
- Electric

Additional Options:

- Stirrups with Adapter
- Monitor Tray Board with Chart Holder
- Push/Pull Bar with *optional* Chart Holder
- Vertical Oxygen Tank Holder

SPECIFICATIONS

Model Number	750-NE
Weight Capacity	750 lbs. (340.2 kg)
Overall Length	88" (223.5 cm)
Overall Width	34" (86.3 cm)
Height Range	
Low	23" (58.4 cm)
High	32.75" (83.2 cm)
Litter Positioning	
Fowler Backrest	0°-90°
Trendelenburg	+12°
Reverse Trendelenburg	-7°
Head Rest Range	
Up	1.5"
Down	1.5"
Total Range	3"
Head Rest Tilt	
Range	+30° - 18°
Patient Surface	82" x 30" (208.2 x 76.2cm)
Side Rails	12.5" x 56" (31.75 x 142.2 cm)
Casters	8" (20.3 cm)
Warranty	1 year parts and labor



✓ **LATEX-FREE**

1

RAIL SYSTEM:

Easy to wipe clean and disinfect in the ongoing battle with germs and bacteria. This 56" anodizing aluminum side rail is infinitely adjustable and can be locked at any height without interfering with the surgical process or patient entry or exit onto stretcher.

2

OXYGEN TANK HOLDER:

Pedigo's ergonomic Quick-Release O₂ Tank Holder is side mounted and conveniently located under the litter to reduce caregiver bending. (US Patent 7,865,984 B2)

3

RETRACTABLE PUSH HANDLES:

Easily extended or retracted push handles by engaging the safety release lock and rotate each handle upward.

4

HYDRAULIC FOOT PEDALS:

Hydraulic Foot Pedals are located on each side of the stretcher within the patient care zone. Pedal function graphics are color-coded for quick and easy recognition.

5

INSTANT STEER:

Mounted at the head end of the stretcher, the unique foot controlled feature places an additional two wheels on the floor for improved traction and six wheel steering.

6

CENTER COLUMN SUPPORT:

Sturdy solid steel column with single hydraulic cylinder design uses sixteen precision ball-bearings to smoothly raise and lower the patient.

